



VIATM Software Theory of Operations

DOCUMENT NUMBER:

CG-00042

DOCUMENT REVISION:

C

Effective Date:

2025-Aug-04

Table of Contents

Revision History	4
Introduction	5
VIA CNV and AOH Segmentation Algorithms	6
Copy Number Segmentation Concepts	6
SNP-FASST2 Algorithm Parameterization	7
BAM MultiScale Reference Concept for NGS	9
OGM Data Integration into VIA	12
OGM Data Workflow for Segmenting CNVs with SNP-FASST3	12
B-Allele Frequency (BAF) from OGM with VIA Software and SNP-FASST3	12
SNP-FASST3 CNV and AOH Performance for OGM	13
VIA CNV and AOH Quality Metrics for OGM Data	15
Concepts of Structural Variants (SV) in VIA	16
Region Types	16
Uniparental Disomy Detection	21
HMM-based Approach to Detect UPD Events	21
Parent of Origin Calculations	23
Error Tolerance	24
Missing Parent Case	24
HRD Genomic Scar Analysis Overview	24
Genomic Instability Scoring for HRD	24
HRD Genomic Scar Processing and Definitions	25
HRD Genomic Scar Performance	26
Aneusomy Detection and Analysis	27

Aneusomy Calling	27
Aneusomy Confidence Scoring	28
Resolving Conflicts Between Whole Chromosomes and Chromosome Arms	28
Expected Workflow for Integrating Aneusomy Calls into the Results Table	30
Guidance for Interpreting Aneusomy Confidence Scores	31
Automatic Pre-Classification Decision Trees	33
Keywords and Syntax For the Decision Tree Rules	33
Basic Functions	35
Classifying an Event	37
Evaluating the type of Event (e.g., copy number gain, AOH, SNV, etc.)	37
Evaluating Event Size and Location	39
Comparing Events to Region lists and Evaluating Similarity Scores	42
Syntax for Structural Variants (SV)	48
Classifying OGM SVs	49
Evaluating Sample Attributes	50
Example Decision Tree Script	51
Tiered Variant Decision Tree for Hematological Malignancies	52
Schematic for Disease-specific Decision Tree	52
Decision Tree Script For Disease Specific Tiering	53
Technical Assistance	59
Legal Notice	62
Patents	62
Trademarks	62

Revision History

REVISION	NOTES
A	Initial document release.
B	Updated for VIA 7.1 software release
C	Updated for VIA 7.2 software release

Introduction

Variant Intelligence Applications™ (VIA) software is a complete and integrated solution for the visualization, interpretation, and reporting of genomic variants from multiple technology types. By supporting multiple genome-wide data modalities, VIA software provides the most comprehensive view of genomic variants of any interpretation, annotation, and reporting software tool available. As a platform-agnostic tertiary analysis solution, VIA stores and manages distinct types of genomic data from various platforms (see **Table 1**) enabling the extraction of meaningful insights from a combined analysis. The software includes algorithms to detect copy number variants (CNV) from major microarray vendors, optical genome mapping (OGM), and next generation sequencing (NGS) methodologies as well as Absence of Heterozygosity (AOH), from data types that assess B-allele frequency. VIA also provides intelligent interpretation assistance to analyze CNVs, Loss of Heterozygosity (LOH) and Structural Variants (SV) from OGM data. As a centralized analysis solution spanning technologies and application areas, VIA software provides an efficient environment to keep pace with advancements in technology while retaining access to historical platform data. By being adaptive to whichever technology is used to generate CNV, LOH, or SV genomic variants, VIA software provides rich annotations for the co-analysis of sequence variants from NGS to provide a complete picture of genomic variation and reveal more answers for disease association.

Table 1. Common platforms supported in the software.

PLATFORM	EXAMPLE ASSAYS	ASSOCIATED FILE TYPES
BIONANO	OGM	ogm.bam
		ogm.vcf
THERMO FISHER/ AFFYMETRIX	Affymetrix arrays output .cel file format	.cel
	CytoScan 750K	.cychp
	CytoScan HD	.cyhd.cychp
	CytoScan XON	.xnchp
	OncoScan	.oschp
	CytoScan HT-CMA, SNP6	.cel
ILLUMINA	CytoSNP12, CytoSNP850K, Infinium Omni, GSA, GSA-Cyto, GDA, GDA-Cyto	.txt (Final Report files)
	Infinium HumanMethylation450, MethylationEPIC	.gtc
	Infinium HumanMethylation450, MethylationEPIC	.idat
AGILENT	SurePrint G3 CGH + SNP Bundle, 4x180K	.txt
	GenetiSure Cyto 4 x 180K CGH+SNP	
NGS	WGS, WES, Panels	CNV = .bam
		Seq Var = .vcf, .vcf.gz
		.json.gz (generated by Nirvana)
CUSTOM	Custom CNV with probe or segment values	.txt (tab delimited)
	Custom Seq Var with annotations	.vcf

The principles and algorithms applied across platforms in VIA software are shared and the flexible parameterization capabilities of the software enable customization to address the nuances of each technology. This document provides a description of the workflows and algorithms fundamental to data processing in VIA with a summary of the software's performance, for example, datasets. OGM data from a VCF or a BAM file can be uploaded into VIA from Access™ or from the VIA homepage. The data is processed using the settings determined in the sample type. This may include application of a decision tree, which assigns pre-classifications to the

variants detected in the sample. Depending on the overlap or similarity with the coordinates of known regions involved in pathogenicity or other regions of interest determined by the user, a decision tree may be applied. Users can visualize and re-classify variants, review links to external databases, write detailed variant interpretation text and export these into report templates. These data are stored in the VIA Server and are accessible for review through the user interface (UI) by other authorized users.

VIA CNV and AOH Segmentation Algorithms

Copy Number Segmentation Concepts

VIA software will arrange the ratios according to their position along the chromosome. Each probe is represented as a small gray dot along the length of a chromosome in the genome and chromosome plots; the user-specified calling thresholds seen as blue and red horizontal lines call certain regions as a **Gain**, **Loss**, **Amplification**, or **Homozygous Loss**. The most trivial calling algorithm would be to simply use these thresholds and the probe locations, but this can cause noisy output, marking any probe that exceeds these limits as a copy number (CN) change event. Different approaches can be utilized, and VIA offers three segmentation algorithms: a circular binary segmentation (CBS)-based, SNP Rank, and two iterations of a hidden Markov model (HMM)-based algorithm, SNP-FASST2 and a new algorithm SNP-FASST3.

SNP-FASST SEGMENTATION

Intensity and BAF information are used in tandem to segment the genome and associate the most likely statehood for the segment. Single Nucleotide Polymorphism Fast Adaptive State Segmentation Technology (SNPFASST) defines states that represent all *possible* combinations of LogR and BAF states in a matrix of ninety-six combinations of copy number and allelic changes.

Effective for detection of mosaic and multi-clonal cases, these probability calculations occur throughout the genome based on the user defined minimum probe floating window. Generally, the higher the number of probes to create a region, the more confidence that the copy number change is real. SNP-FASST is based on the Hidden Markov Model, but the defined states are Normal, Gain, Amplification, Loss, and Homozygous Loss as seen in **Figure 1**. The state boundaries are adaptive, meaning the user can adjust the thresholds (blue/red lines) so that the algorithm adheres to the data.

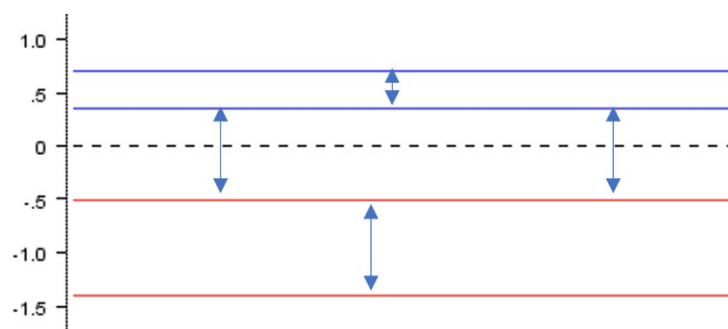


Figure 1. CN state threshold setting for SNPFASST region calling.

SNP-FASST2 Algorithm Parameterization

Parameterization of the segmentation enables flexibility to achieve desired sensitivity and adapt to technology nuances, which are set within the VIA administer panel. The key parameters impacting segmentation for copy number events are:

- Copy number state thresholds for each CN event type, and corresponding settings for autosomes and sex chromosomes wherein the tighter thresholds will increase sensitivity
- Significance Threshold is a critical value for SNP-FASST2 segmentation as functionality and influence on region calling. Increasing the stringency of the p-value to segment the genome will enhance CNV and AOH calling from data patterns with higher confidence.
- The minimum number of probes to make a CN segment
- Maximum Contiguous Probe Spacing (in Kilo base-pairs - Kbp) is a parameter used to limit the segments such that if there are any two neighboring probes that are separated from each other by more than the specified distance, the segmentation algorithm will stop at the last probe location. This will result in no calls being made in these areas. A good example would be the centromere loci.

This algorithm identifies trends in the B-Allele Frequency for copy neutral events and is highly informative in identifying long contiguous stretches of homozygosity (LCSH), iso uniparental disomy (iUPD), chimerism, or maternal cell contamination (MCC). **Figure 2** illustrates these unique features as described below.

- **Minimum LOH length** - defines the size of the contiguous AOH; it is a key SNP calling parameter.
- **Homozygous Threshold** – defines a homozygous state for the SNPs between the yellow and black lines.
- **Heterozygous Threshold** – defines an SNP between the purple and yellow lines and adjusts sensitivity for calling AI.

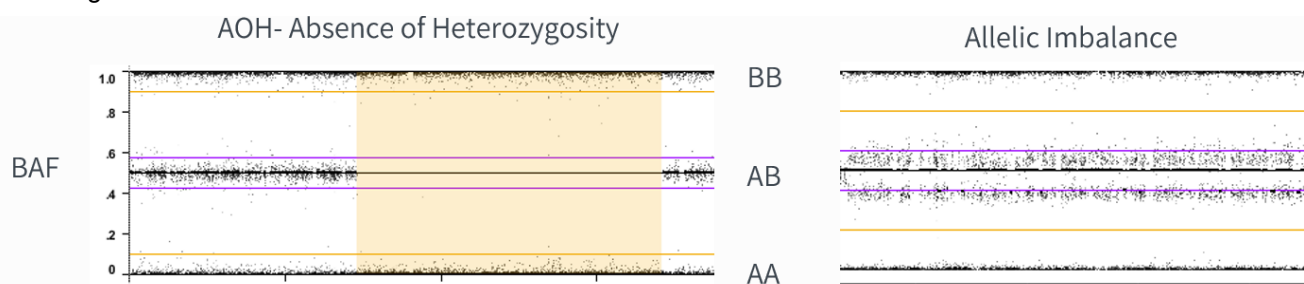


Figure 2. SNP- Allelic Imbalance/AOH Detection

SNP-FASST3 SEGMENTATION

SNP-FASST3 segmentation algorithms build on the previous version, SNP FASST2, and offers improved segmentation particularly for mosaic events. The main enhancements are:

- Support for mosaic calling
- Improved calling of events when probes fall on the one-copy loss or one-copy gain threshold lines
- Integration of BAF and Log R for copy number calling with SNP arrays (for SNP-FASST3)

This is achieved through the addition of mosaic states with the processing settings and a modified probability curve for optimal region segmentation. The distribution used, Plateau Pseudo Distribution (PPD), is a modified normal distribution with a plateau in the center, as seen in **Figure 3**.

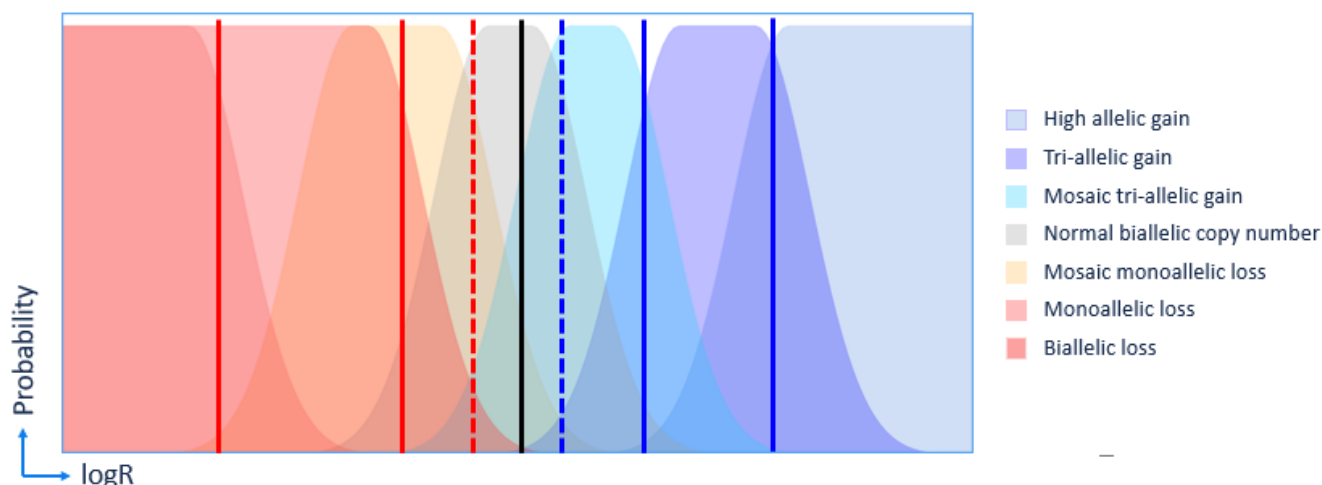


Figure 3. An example Plateau Pseudo Distribution (PPD) plot for loss states.

The addition of mosaic states to the Log R and BAF states, which represent a state between normal and adjacent states allows use of a single settings profile to call copy gains and losses as well as mosaic gains and losses. Entering these mosaic states requires orders of magnitude more significant than entering non-mosaic states (e.g., if the significance threshold for non-mosaic is set to 1E-8 then, by default, the threshold for the mosaic state would be 1E-20. This offset can be modified by the user in the **Processing** settings).

The algorithm also fills gaps between the loss states and the gain states such that the probability distribution function for homozygous loss/big loss, loss, and mosaic loss (and similarly for high gain, amplification/gain, and mosaic gain) will not have dips between the distributions. This improves calling when probes fall along the threshold lines where these dips were previously located.

In addition, SNP-FASST3 incorporates Log R and BAF data together for single-copy loss and mosaic copy loss states to improve calling, where previously only Log R values were used to determine copy number. This makes the algorithm less likely, for instance, to make a single-copy loss call in a region where the BAF track looks heterozygous. Values for gain and loss are 1/6 of the one-copy gain/loss thresholds respectively and for mosaic allelic imbalance, it is $(2.5 + \text{allelic imbalance threshold})/6$.

Several parameters that affect the behavior of mosaic calling can be adjusted by the admin in **Processing** settings for SNP-FASST3 algorithm:

- Minimum mosaic threshold (%) – level of mosaicism to be detected, corresponding to the desired aberrant cell fraction
- Mosaic CN significance offset – amount by which to increase stringency of the significance threshold for mosaic CN states
- Mosaic SNP significance offset – amount by which to increase stringency of the significance threshold for the mosaic imbalance state

SNP-FASST3 also decouples the significance threshold, which is the most impactful processing setting, for CN and AOH segmentation to allow for additional parametrization to achieve desired sensitivity for segmenting CN and AOH independently.

RANK AND SNP RANK SEGMENTATION ALGORITHM

The Rank Segmentation algorithm is a statistically based algorithm, similar in concept to the Circular Binary Segmentation (CBS) algorithm developed by Adam Olshen at Sloan-Kettering Institute (Olshen AB, Venkatraman ES, Lucito R, Wigler M. Circular binary segmentation for the analysis of array-based DNA copy number data. *Biostatistics*. 2004 Oct; 5(4):557-72). The CBS algorithm has been modified to significantly improve processing speed by using a normal distribution function to assess for change points as opposed to the non-parametric permutation-based statistics used in the original CBS algorithm. The result is segmentation of the genome into clusters of uniform ratios. A recursive algorithm, the genome is continuously divided into smaller and smaller units until no region can be further segmented. Significance Threshold is the single parameter that controls whether a region is to be segmented out or not. The logic is to start by rank ordering the log-ratio probe values, consider the distribution of the probe ranks in a region and then compare this information to the distribution of probes in the adjacent segment (to the left and right). If these distributions are significantly different, generating a significance value less than the Significance Threshold, then the segments are divided. The process stops if no segment can be found in an interval that is significantly different than its neighbors. At completion, the entire genome can be represented as a series of segments, each having a cluster value which is the median log-ratio value of all the probes in that region, plotted as horizontal black lines. The calling algorithm then uses the cluster values and the user-defined thresholds to establish regions of copy number variations. The SNP Rank Segmentation is like the Rank Segmentation algorithm but also considers B-Allele frequency values from SNP arrays to segment the genome. Using the B-allele frequency in conjunction with log ratios (providing copy number results) allows for better segmentation.

For the Rank, SNP Rank, FASST2, and SNP-FASST2 Segmentation algorithms, a significance threshold needs to be set in the Analysis panel so the sensitivity can be adjusted. The smaller the number, the less sensitive the algorithm is in creating a new segment. So, if some known aberrations are not being called because they are too small, this value should be increased. This setting is inversely proportional to the number of probes: the larger the number of probes, the smaller the value used for this setting, ensuring valid results. Many probes at a setting of 1E-6 or lower have been processed.

BAM MultiScale Reference Concept for NGS

The BAM MultiScale Reference (MSR) method functions well with both shallow and targeted sequencing data as well as WGS/WES with normal depth of coverage. It builds a reference file from a set of normal samples and uses adjustable dynamic binning. The method uses a Hidden Markov Model to segment the genome into target areas using the reads in targeted regions and the backbone areas using the off target reads and additional areas. Coarse binning is used in the backbone areas to provide the copy number baseline as well as large copy number events and fine binning is used in target areas to provide high resolution copy number detection in targeted regions. The adjustable dynamic binning is very flexible allowing adjustment of the minimum bin width based on the depth of coverage. The dynamic binning allows the target regions to get more coverage and the backbone regions, less coverage but the backbone still gets coverage. **Figure 4** illustrates the dynamic binning approach.

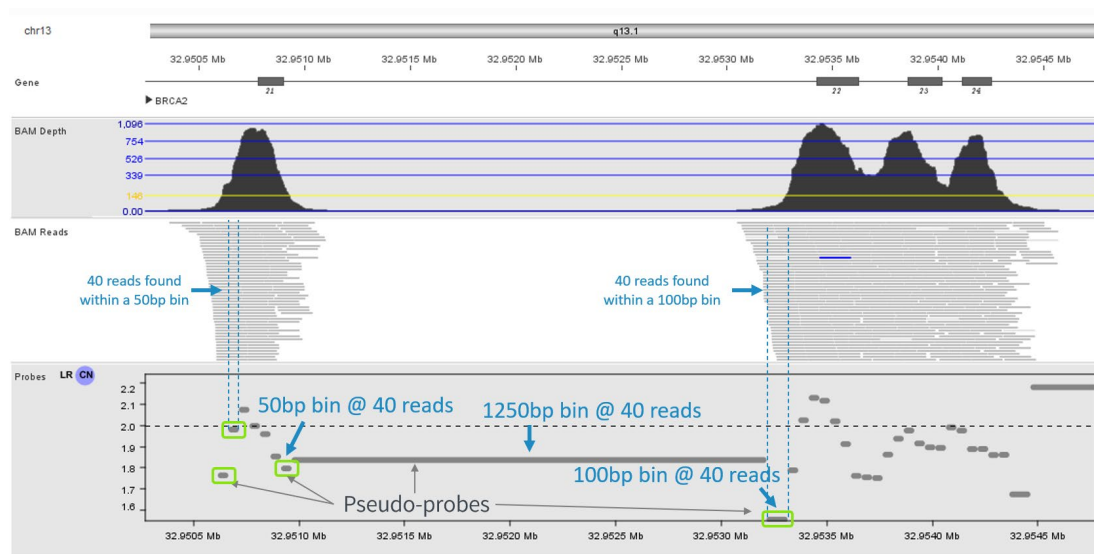


Figure 4. Bin formation using the BAM MSR Builder tool

The BAM MSR method is a read-depth method that uses a pooled reference file to generate pseudo-log ratios based on the reads that it packages as a separate utility, the BAM MultiScale Reference Builder program. It also generates B-allele frequencies based on the reads at SNP locations. After logR bins are generated from the NGS methods, segmentation for CNV and AOH/LOH is performed using any of the segmentation algorithms available in VIA. Instructions on the parameterization of the BAM MSR is provided in the *VIA Software User Guide* (CG-00043).

SNP-FASST3 CNV PERFORMANCE FOR WGS (WHOLE GENOME SEQUENCING)

Copy number variant detection performance of SNP-FASST3 using whole genome sequencing data was assessed for both MultiScale Reference (MSR) and self-reference calling methods using simulated data. For the MSR analysis, the applied MSR file was established with the following settings: minimum bin width = 50, average read length = 50, target reads per bin = 300. Copy number gains and losses were simulated at twelve size ranges between 5 kb to 3.5 Mb at target allele fractions of 10%, 20%, 30%, 40% and 50%. Simulated variants were incorporated into a baseline Illumina WGS sample sequenced to a depth of 45X. Residual CNV calls in the actual sample were identified with SNP-FASST3 and then subtracted from the final results so that performance evaluation was done using only simulated events. Five gains and five losses were simulated for each size bin and variant allele fraction for a total of 600 events with no more than ten events in each simulated sample. Resultant BAM files were analyzed with SNP-FASST3 and CNV calls were assessed for calling accuracy. Centromeric and annotated segmental duplication regions were excluded from the analysis. CNV calls were made using the self-reference method for coverage normalization as well as with an MSR composed of six WGS samples sequenced to a comparable depth with the same methods. See **Figure 5** and **Figure 6**.

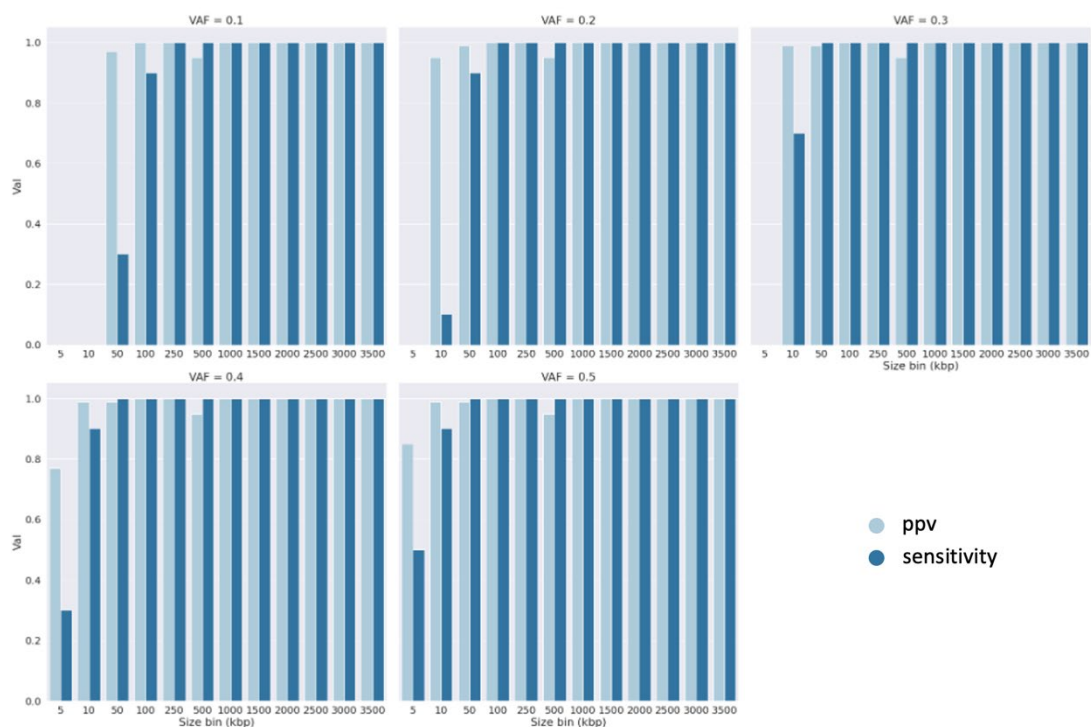


Figure 5. Sensitivity and PPV of simulated CNV events detected using SNP-FASST3 self-reference method

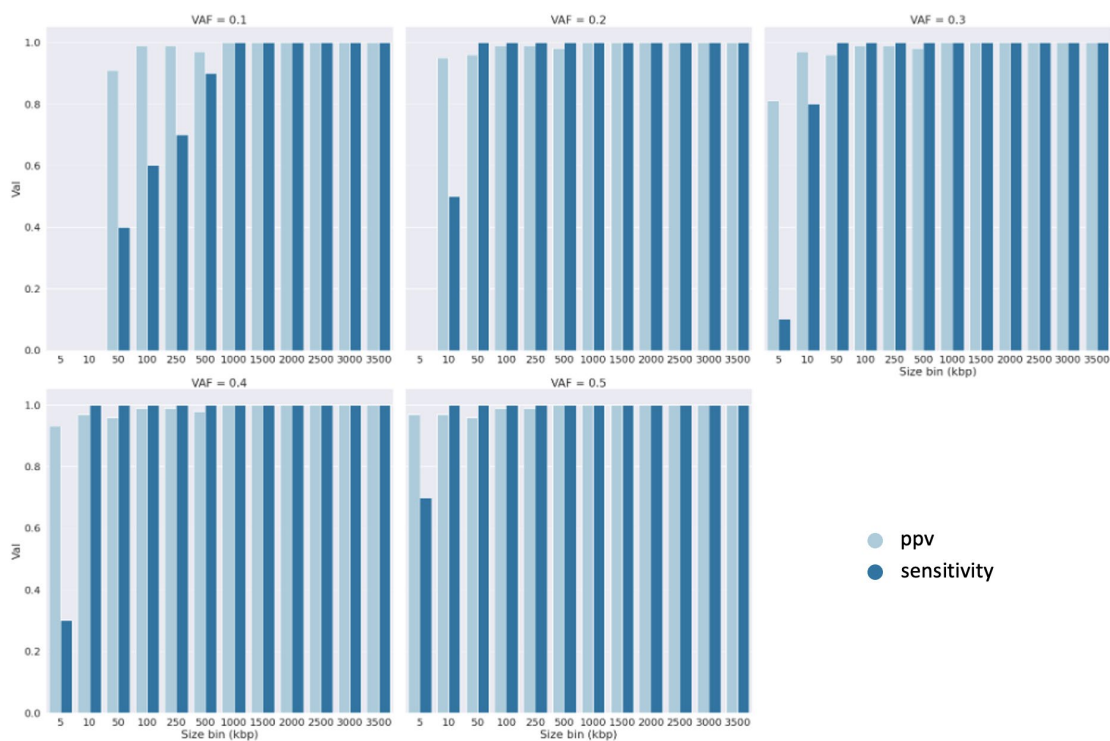


Figure 6. Sensitivity and PPV of simulated CNV events detected using SNP-FASST3 MSR method

OGM Data Integration into VIA

OGM Data Workflow for Segmenting CNVs with SNP-FASST3

OGM samples can be uploaded directly from Access by selecting the option to 'Upload to VIA' when submitting a sample for processing or the data can be manually uploaded into VIA through the Data or the Batch Import method leveraging the ogm.bam (with the accompanying ogm.bam.bai) and ogm.vcf files for each sample.

It is recommended to process OGM BAM files with an OGM BAM MultiScale Reference matching the effective coverage, sex, and genome build. The OGM MSR files are created with a set of cytogenetically normal samples through the BAM MultiScale Reference Builder. Like the BAM MultiScale method for NGS, a set of normal OGM BAM data files are used to construct dynamic bins of the coverage profile for the reference dataset to generate a MSR file inclusive of bin positions with expected coverage levels. The MSR file is applied as an *in silico* reference for an experimental data file to measure coverage difference at each bin for CNV segmentation processing. Bionano has generated shareable OGM MSR reference files from male and females control samples at 80x, 160x, and 300x effective coverage levels for each genome build leveraging the settings seen in **Table 2**.

Table 2. OGM MSR Reference Builder Settings

OGM Effective Coverage	80x	160x	300x
Minimum Bin Width (bp)	1000	1000	1000
Average Read Length (bp)	100000	100000	100000
Target Reads per Bin	80	200	200
Maximum Neighbor Bin Gap (bp)	100	100	100

B-Allele Frequency (BAF) from OGM with VIA Software and SNP-FASST3

The B-allele frequency of a label is calculated as the ratio of OGM molecules observed with label fluorescence to the total number of molecules aligned to the position. Heterozygous labels are defined as label motif sites that overlap heterozygous SNPs, where a fraction of molecules is observed with a lack of fluorescence at that point. Most labels are homozygous-present (each label overlaps with six homozygous normal base pairs). There is a 3-step process for identifying a BAF data point:

1. The data gets filtered to labels that overlap a known SNP (MAF > 5%).
2. Labels are further filtered for those that cluster into three well-separated clusters (hom-absent, het, hom-present) when observing across 180 control samples, based on the Silhouette score (**Figure 7**).
3. BAF values are normalized based on BAF distributions in control samples.

BAF data point generation and AOH calling performance has been established with OGM datasets at 300x effective coverage. Performance for calling aberrations based on BAF data, such as regions of AOH and polyploid genomes, has not been verified on lower coverage datasets. Please see *CG-30548-Visualizing-Different-Classes-of-Structural-Variants-in-Bionano-Access-Software.pdf* for guidance on triploidy assessment in Bionano Access for constitutional samples.

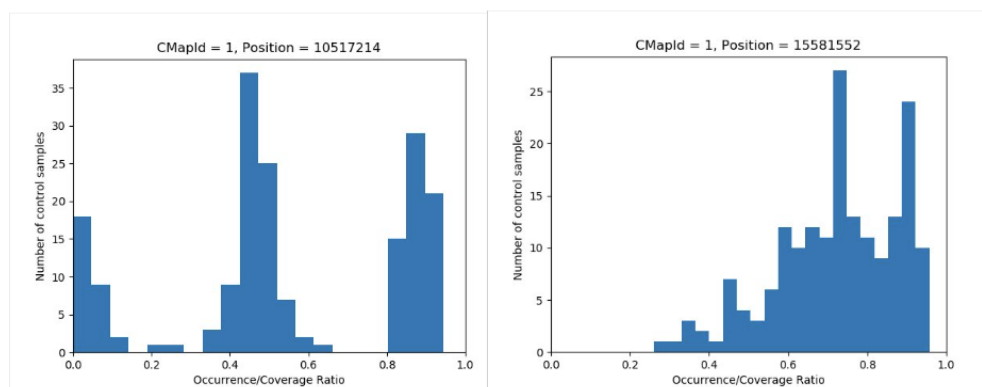


Figure 7: Example clusters for a good (left) and bad (right) labels with a SNP loci

Event significance scores for CNVs

The event significance is intended as a tool for filtering out CNV calls that are likely to be false positives, based on the size of the event and the estimated allele fraction/copy number. Statistically, it is the probability that the event is a false positive given that the only source of noise in the region is background noise due to random molecule sampling. Noise that occurs at a much larger scale (related to DNA packing, for example) or noise in regions of high variation such as centromeric and telomeric regions is not modeled in the event significance score. As a result, events that are called due to these more systematic noise sources will have unreliable values for event significance and are likely to be reported as highly confident due to their large size and strong signal, even though they are likely to be false positives. Therefore, the event significance will only reflect the overall event significance in clean samples (without artifacts from DNA packing) and in regions that are typically reliable (e.g., not telomeres and centromeres). An event score threshold of $1e-8$ was applied for the determination of CNV performance for OGM data.

SNP-FASST3 CNV and AOH Performance for OGM

CNV and AOH calling performance was evaluated in samples with known events, as well as simulated samples. For AOH evaluation, 230 AOH/LOH events of sizes ranging from 1Mb to 100Mb were simulated on 300x datasets at various aberrant cell fractions as low as 5% aberrant cell fraction (ACF), for a total of 1,350 AOH/LOH events evaluated. The processing settings used for AOH evaluation are specified in the Appendix. Sensitivity and precision were calculated for each size range and cell fraction separately. It was observed that recall for 20-25 Mbp AOH events is 92% at 25% ACF. Additional validation was performed using a cohort of constitutional and cancer samples for which orthogonal testing had been performed and for which cell counts and LOH events were available. In 15 samples containing 37 known AOH events, all events were called in 14/15 samples. One 14.6 Mbp AOH event was not called; other AOH events that were greater in size were called. Most false positive calls were under 10 Mbps, and false positive calls in the 40-10 Mbps range could be distinguished using manual review. See **Figure 8**.

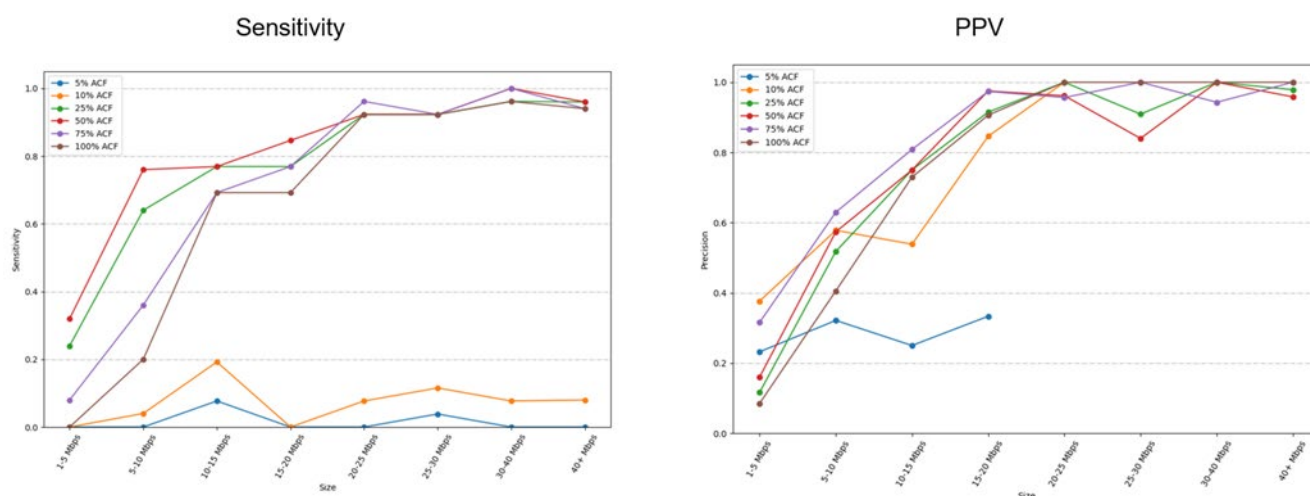


Figure 8. AOH calling performance using simulated data at 300x effective coverage. True positives are defined as 50% overlap required between a simulated event and an AOH call.

For CNV evaluation evaluated using default processing settings, 280 CNV events ranging in size from 175kb to 8Mb were simulated on the 22 autosomal chromosomes at multiple coverage depths and allele fractions. Data, seen in **Figure 9**, **Figure 10**, and **Figure 11**, were simulated to represent 50% VAF at the 400 Gbp (80x) level, 20%, 30% and 50% at the 800 Gbp (160x) level and 5%, 10%, 20%, 30% and 50% at the 1.5 Tb (300x) level. Sensitivity and PPV were calculated for each size range and cell fraction separately. True positives were defined as events with 80% overlap required between simulated event and CNV call. Variants overlapping the OGM CNV Mask region by 45% or more were excluded from analysis. For a detailed description of the creation and content of the CNV mask, see *Bionano Solve Theory of Operation: Structural Variant Calling* (CG-30110). For a summary of CNV performance, see *Bionano System Application Specifications* (CG-00008).

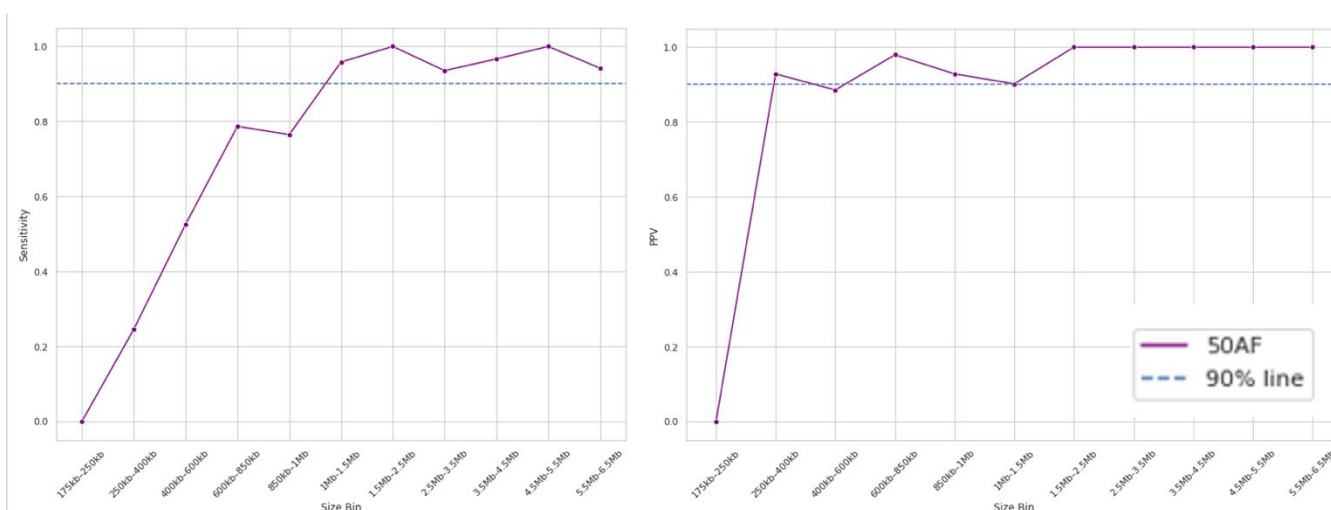


Figure 9. CNV detection performance at 400 Gbp/80x coverage

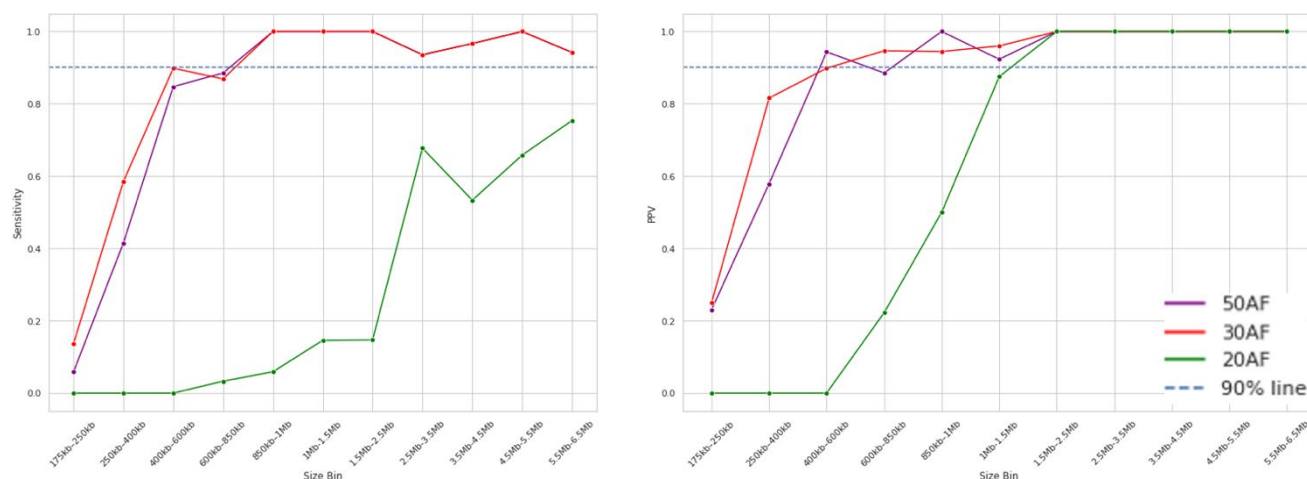


Figure 10. CNV detection performance at 800 Gbp/160x coverage

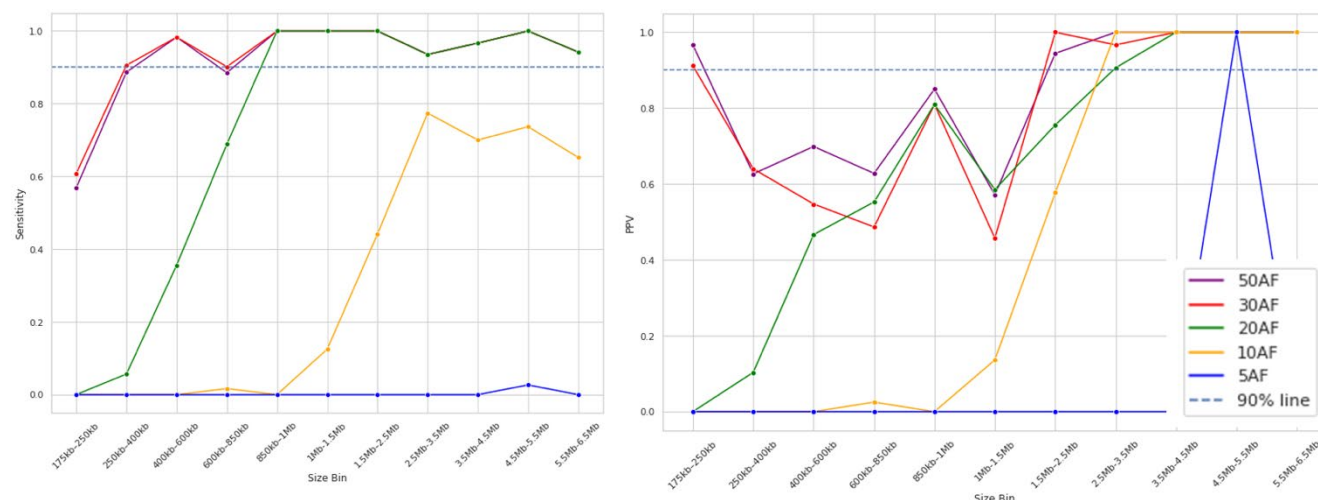


Figure 11. CNV detection performance at 1.5 TB/300x coverage

VIA CNV and AOH Quality Metrics for OGM Data

VIA includes software-specific metrics for assessing the data quality to inform reliability in the CNV and AOH called conducted by the software algorithms. These metrics are described below with typical ranges for good quality OGM data from normal samples. Note that the quality metrics can be influenced by samples with complex and numerous chromosomal aberrations, which should be factored into consideration when reviewing a sample's quality metrics.

CN Quality: A score representing the probe-to-probe variance measuring on average how much successive probes differ from each other. It can help indicate reliability of the data; a higher CN Quality score indicates higher variance and less reliable data while a lower CN score indicates cleaner, more reliable data. A high score may indicate DNA quality or other issues during wet lab prep.

- Typical CN Quality score for good quality OGM data (300x): <0.005
- Typical CN Quality score for borderline (warning) quality OGM data (300x): >=0.005

- Typical CN Quality score for poor quality (failing) OGM data(300x): ≥ 0.015

BAF Quality: Represented as the percentage of SNP probes with BAF values in the allelic imbalance region of the BAF track. This quality score indicates how well data points on the BAF plot are clustering. It reports the percentage of SNP probes with BAF values between the heterozygous imbalance threshold and homozygous value threshold as defined in the **Processing** settings.

- Typical BAF Quality score for good quality OGM data: < 0.01
- Typical BAF Quality score for borderline (warning) quality OGM data (300x): ≥ 0.01
- Typical BAF Quality score for poor quality (failing) OGM data(300x): ≥ 0.1


Concepts of Structural Variants (SV) in VIA

Region Types

SV Events in VIA are typically defined by multiple breakend positions plus confidence region with relation to the reference sequence. The region types are combinations of these positions/regions as illustrated in **Figure 12**.

FULL EXTENT REGION

This region is defined by the most upstream point of the confidence region of the start position, and the most downstream point of the confidence region of the end position. This is typically used for Dosage Effect events.

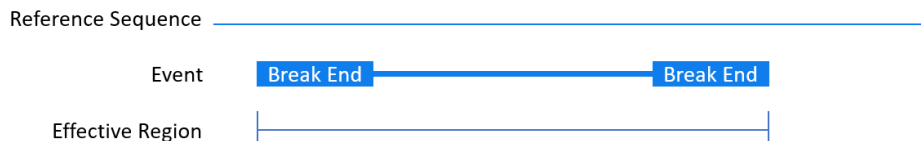

BREAK END REGIONS

This type of region is defined by only the region of the break ends and their confidence interval. This is typically used for gene disruption events, as well as translocation and fusions.



REGION PADDING

Padding may be added to the regions defined above as needed.

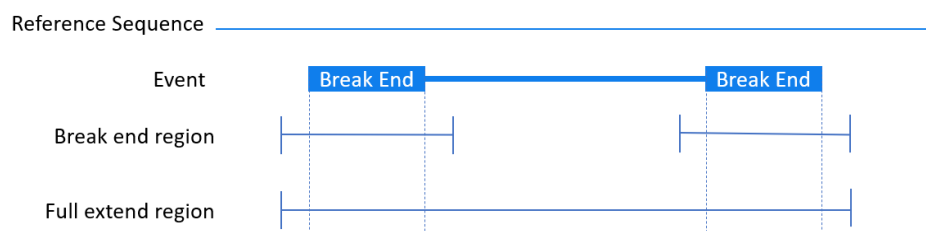


Figure 12. Regions

SV EVENT TYPES

Definitions

- Visualization Region: Region used to display the SV event in the SV Track in VIA.
- Column Region: Region displayed in the sample event table, variant details tab. This region is also used for region overlap calculations for various annotation regions/lists in the sample event table.
- Variant Details Gene Table: Region used to generate the gene table in the variant details tab.

DUPLICATION

This is considered a Dosage Effect event.

- Visualization: Full Extent Region
- Column Region: Full Extent Region
- Gene Panel Evaluation: Full Extent Region
- Variant Details Gene Table: Full Extent Region

DELETION

This event is considered a Dosage Effect event.

- Visualization: **Full Extent** Region
- Column Region: **Full Extent** Region
- Gene Panel Evaluation: **Full Extent** Region
- Variant Details Gene Table: **Full Extent** Region

INSERTION

This event is considered a Dosage Effect event.

- Visualization: **Full Extent** Region
- Column Region: **Full Extent** Region
- Gene Panel Evaluation: **Full Extent** Region
- Variant Details Gene Table: **Full Extent** Region

INVERSION

This event is considered a Gene Disruption event.

- Visualization: **Full Extent** Region
- Column Region: **Break end** Regions
- Gene Panel Evaluation: **Break end** Regions
- Variant Details Gene Table: **Break end Regions** plus padding

INTRACHR FUSION

This event is considered a Gene Disruption event.

- Visualization: **Break end** Regions
- Column Region: **Break end** Regions
- Gene Panel Evaluation: **Break end** Regions
- Variant Details Gene Table: **Break end** Regions plus padding

INTERCHR TRANSLOCATION

This event is considered a Gene Disruption event.

- Visualization: **Break end** Regions
- Column Region: **Break end** Regions
- Gene Panel Evaluation: **Break end** Regions
- Variant Details Gene Table: **Break end** Regions plus padding

INVERTED DUPLICATION

This event is considered both Dosage Effect and Gene Disruption event.

- Visualization: **Full Extent** Region
- Column Region: **Full Extent** Region
- Gene Panel Evaluation: **Full Extent** Region
- Variant Details Gene Table: **Full Extent** Region plus padding

SV LENGTH

For all SV types other than Inversion, the value displayed in the length column in VIA is the absolute value of the SV length reported in the OGM.VCF file. In the case of Inversion, it is the size of the inversion (SVLEN in the OGM VCF would be 0 because there was no change in actual size), which is calculated using the POS and End specified in the OGM VCF file.

ISCN FORMATTING FOR SV IN VIA

Table 3 and **Table 4** summarize the positional placement of coordinates used with the nomenclature according to the ISCN for OGM data types.

Table 3. ISCN intended coordinates

Event Name	BNGTYPE	Description of Start and End
Deletion	deletion	First and last base of deleted sequence
Tandem Duplication	duplication	First and last base of duplicated sequence
Inverted Duplication	duplication_inverted	First and last base of duplicated sequence
Insertion	insertion	Base before and base after inserted sequence
Inversion	inversion_paired	First and last base of inverted sequence
Inversion Breakpoint	inversion_partial	First and last base of inverted sequence
Interchr Translocation	translocation_interchr	Position of each break-end
Intrachr Fusion	intrachr_fusion	Position of each break-end

Table 4. ISCN coordinates relative to the VCF

Event Name	BNGTYPE	ALT	Start	End
Deletion	deletion		POS + POS.Cl.min + 1	POS + POS.Cl.max
Tandem Duplication	duplication	<DUP:TANDEM>	POS + POS.Cl.max + 1	END+ END.Cl.min
Inverted Duplication	duplication_inverted	<INV>	END + END.Cl.min	POS + POS.Cl.max
Insertion	insertion	<INS>	POS + POS.Cl.min	POS + POS.Cl.max
Inversion	inversion_paired	<INV>	POS + POS.Cl.min + 1 ~ POS + POS.Cl.max + 1	END + END.Cl.min ~ END + END.Cl.max
Inversion Breakpoint	inversion_partial	N].]	POS + POS.Cl.min + 1 ~ POS + POS.Cl.max + 1	Mate.POS + MATE.POS.Cl.min ~ Mate.POS + MATE.POS.Cl.max
Inversion Breakpoint	inversion_partial	[.[N	POS + POS.Cl.min ~ POS + POS.Cl.max	Mate.POS + MATE.POS.Cl.min - 1 ~ Mate.POS + MATE.POS.Cl.max - 1
Interchr Translocation	translocation_interchr	N].] or N[.[or].]N or].]N	POS + POS.Cl.min ~ POS + POS.Cl.max	Mate.POS + Mate.POS.Cl.min ~ Mate.POS + Mate.POS.Cl.max
Intrachr Fusion	intrachr_fusion	N].] or N[.[or].]N or].]N	POS + POS.Cl.min ~ POS + POS.Cl.max	Mate.POS + Mate.POS.Cl.min ~ Mate.POS + Mate.POS.Cl.max

TABLE KEY

- CI = confidence interval for primary breakpoint
- min = left or start position of confidence interval (this will be a negative value)
- max = right or end position of confidence interval
- POS = primary breakend position
- END = secondary breakend position
- Mate = linked or partner breakpoint

NOTES

- **NOTE:** The current convention is to list the end before the start for inverted duplications.
- The VCF orders SVs by position
- The ALT field for Inversion Breakpoint, Interchr Translocation and IntraChr Fusion indicates the orientation of the breakends as specified in the VCF v4.2 format specifications. More information on VCF formatting is available in the *OGM File Format Specification Sheet* (CG-00045).

ISCN EXAMPLES

Table 5 provides example nomenclature for different types of SVs.

Table 5. Example nomenclature for different types of SVs

Event Name	BNGTYPE	ISCN
Deletion	deletion	ogm[GRCh38] 1p21.1(103,607,985_103,682,039)x1
Tandem Duplication	duplication	ogm[GRCh38] dup(7)(q11.21q11.21)(62,926,148_62,984,222)
Inverted Duplication	duplication_inverted	ogm[GRCh38] dup(15)(q13.1q13.1)(28,486,912_28,327,161)
Insertion	insertion	ogm[GRCh38] ins(1;?)(p21. 1;?)(103,935,603~103,948,580;?)
Inversion	inversion_paired	ogm[GRCh38] inv(1)(p21.1p21.1)(103,607,985~103,703,390_103,764,118~103,772,456)
Inversion Breakpoint	inversion_partial	ogm[GRCh38] inv(1)(p21.1p21.1)(103,607,985~103,703,390_103,764,118~103,772,456)
Interchr Translocation	translocation_interchr	ogm[GRCh38] t(1;7)(q32.1;p21.3)(202,731,409~202,755,096;7,757,501~7,765,040)
Intrachr Fusion	intrachr_fusion	ogm[GRCh38] 20::20(q11.21::q13.2)(32,497,427~32,502,975::53,552,309~53,553,376)
Loss	loss	ogm[GRCh38] 1p36.33(710,374_711,817)x0~1
Gain	gain	ogm[GRCh38] 1p36.33(710,374_711,817)x2~3

Uniparental Disomy Detection

VIA software allows for the automated detection and analysis of uniparental disomy (UPD) events from either SNP arrays or NGS data. The identification of both isoUPD and heteroUPD are important for genetic interrogations. From SNP arrays and NGS data with corresponding SNP data, VIA software enables laboratories to detect the presence of both isoUPD and heteroUPD when a parental sample is available.

VIA software compares the SNPs of the proband against its parents using an HMM to identify regions where uniparental inheritance is more likely than usual biparental inheritance. For example, SNP positions where the proband is homozygous opposite from one of the parents makes uniparental inheritance, from the other parent, more likely. Regions of isoUPD and heteroUPD are flagged on the proband sample to enable the reviewer to make an interpretation of genes impacted by the aberrant parental inheritance pattern.

HMM-based Approach to Detect UPD Events

EMISSION PROBABILITIES:

We must first model the emission probabilities of the set of genotypes G we have observed given the various states S (hUPD-mom, hUPD-dad, isoUPD-mom, isoUPD-dad, normal). The likelihood of observing a set of genotypes given a particular state S can be defined as the following product over all n positions:

$$P(G|S) = \prod_{i=0}^n P(g_i|S)$$

NOTE: In the log space, we can take the above sum instead. The brunt of the approach is in modeling the individual likelihoods $P(g_i | S)$. Note that we use the general term g which includes all three genotypes at the position - mom (m), dad (d), and proband (p). Remembering that $P(g_p, g_m, g_d) = P(g_p | g_m, g_d) P(g_m, g_d)$ using the chain rule, we can first define $P(g_p | g_m, g_d)$:

Take, for example, $P(AA | AA, BB)$. This quantity will be different depending on what inheritance model we assume. For example, under an assumption of normal Mendelian inheritance, this genotype is impossible. It is also impossible under either of the two paternal UPD models (hUPD-dad and isoUPD-dad). However, under both maternal UPD models (hUPD-mom and isoUPD-mom), $P(AA | AA, BB) = 1.0$. This quantity can be similarly estimated for all possible combinations of genotypes given our theoretical understanding of the five possible inheritance models.

Next we must define $P(g_m, g_d)$. In this case, we empirically estimate from data the probability of all possible mother-father pairs of genotypes in a platform-specific manner.

TRANSITION PROBABILITIES:

We compute transition probabilities based on a theoretical empirical dataset generated from the following four chromosomal archetypes:

- Unaffected - no events
- Affected - single whole-chromosome UPD event
- Affected - segmental UPD alternating between same-parent types, with 10 segments per chromosome
- Affected - lone partial UPD event

The total proportion of affected chromosomes (a) and unaffected chromosomes (u) is computed based on the user-defined significance threshold according to the following formula:

$$a = 2^{\log_{10} t/S}$$

$$u = 1 - a$$

where t is the user-defined significance threshold, and S is a constant currently hardcoded to $1e-5$. It represents the threshold corresponding to 100% affected chromosomes. However, we cap the affected chromosome result at 0.5, to ensure some minimum representation of unaffected chromosomes even with a very lax significance threshold. From this theoretical dataset of chromosomes, we can compute the probability of transitioning at a given position from one of the five inheritance states to any of the others.

Parent of Origin Calculations

The concept of an "informative" probe, which is an SNP probe that is only consistent if inherited from one parent and not the other (for example, when the proband is AA and the mother is AA and father is BB) can be used to effectively determine parent of origin for trio analyses. This means that without any Mendelian error, the probe must have been inherited from the mother. The parent of origin call was made based on which parent a significantly higher number of informative probes favoring inheritance from them. While this works well for trios, for duos, if the Parent of Origin is the present parent, no probes will be informative without knowing the genotype of the missing parent. Moreover, if there is even a single erroneous probe indicating the missing parent, the event may be called (in error) for the missing parent.

Therefore, VIA software implemented a statistically more complete solution to calculate the likelihood ratio for the Parent of Origin of an event using all the probes in the event. To do this, we take each probe in the event and compare the genotype of the proband against the genotypes of each parent and calculate probability that the proband genotype was inherited from that parent. This happens for all the probes in the event and a summary statistic is calculated that indicates how many times more likely it is that the event was inherited from one parent versus the other. The likelihood ratio threshold for calling an event from one parent or another is 10x. This means that for the software to assign a Parent of Origin, the probability of the event having been inherited from that parent must be at least 10 times more than inheritance from the other parent. Hence, the likelihood ratio using all the probes in the event provides a better determination of the parent of origin, especially for events that may not have many "informative" probes, as in the case for duos.

Computing the likelihood of an event origin:

The posterior probability of the event state S (whether it was inherited from the mother or father) given the set of trio genotypes G we are observing can be expressed using Bayes law:

$$P(S|G) = \frac{P(G|S) * P(S)}{P(G)}$$

To avoid evaluating the general probability of observing all trio genotypes in the event, $P(G)$, the Bayes factor is used to determine the marginal likelihood of our two competing hypotheses (mom origin vs dad origin). In this case, we must only compute the likelihood ratio K :

$$K = \frac{P(G|S_m)}{P(G|S_d)}$$

If K is large (over 100, for example), then we can confidently assume maternal origin. If K is small (under 0.01), we can confidently assume paternal origin.

To compute K , only the probability of a set of observed genotypes G given the two possible origin states S (mom, dad) is used. The likelihood of observing a set of genotypes given a particular state S can be defined as the following product over all n positions:

$$P(G|S) = \prod_{i=0}^n P(g_i|S)$$

Note that in the log space, we can take the above sum instead. The brunt of the approach is in modeling the individual likelihoods $P(g_i|S)$. Note that we use the general term g which includes all three genotypes at the position - mom (m), dad (d), and proband (p). Remembering that $P(g_p, g_m, g_d) = P(g_p|g_m, g_d) P(g_m, g_d)$ using the chain rule, we can first define $P(g_p|g_m, g_d)$. This quantity is straightforward to compute given an assumption of normal Mendelian inheritance. For example, we know that $P(AB|AB, AA) = 0.5$. Similar theoretical values can be obtained for all other possible combinations of genotypes. $P(g_m, g_d)$, on the other hand, must be estimated empirically over all possible mother-father pairs of genotypes, in a platform-specific manner.

Error Tolerance

After obtaining the joint probability $P(g_p, g_m, g_d)$ for each possible state, error tolerance is added to each probability due to measurement error. Assuming an error weight of $E = 0.03$, the uniform probability over all 27 possible genotypes is defined as $p_{\text{uniform}} = 1 / 27$. The error-adding function is as follows:

$$p_{\text{new}} = (1 - E) * p_{\text{old}} + E * p_{\text{uniform}}$$

Missing Parent Case

For duo samples with a missing parent, rather than taking the simple product to obtain the joint probability: $P(g_p, g_m, g_d) = P(g_p|g_m, g_d) P(g_m, g_d)$, the sum over all possible genotypes is taken for the missing parent, where a is the available parent, like so:

$$p(g_p, g_m, g_d) = \sum_{g=AA,AB,BB} p(g_p|g_a, g) p(g_a, g)$$

HRD Genomic Scar Analysis Overview

Genomic Instability Scoring for HRD

Homologous recombination deficiency (HRD) is the inability to repair double-stranded DNA breaks using the Homologous Recombination Repair (HRR) cellular pathway, which consequentially results in an acquired chromosomal breakage. Clinical research has shown that cells with HRD are more sensitive to certain therapies and a measurement of HRD can be an effective pharmacogenetic biomarker across various tumor types. To provide a functional evaluation of HR status, HRD genomic scarring is an analysis approach to assess three specific quantifiable signatures of HRD genomic instability. VIA includes a measurement of these three genomic scars to aid with HRD status assessment in cancer samples across technology types.

HRD Genomic Scar Processing and Definitions

Within the Admin section for Sample Types that are set to an Oncology Test Type, selecting the automated **Perform Genomic Scar Calculation** checkbox will activate the analysis during the processing for all associated samples. In brief, genomic CNV and AOH profiles generated across data types are analyzed for scar characteristics through the implementation of three processing steps:

MERGING

During this step, copy number and allelic events are converted into a single genome representation. The two tracks are merged, preserving all existing breakpoints, and the result is stored with no information loss. Each unique combination of CN and allelic calls gets its own merged event state in the resulting single genome representation.

SMOOTHING

The resulting merged track combines similar event types across small gaps as well as the centromere. Each of the three scores has a custom smoothing procedure, which considers each individual score's preferences for minimum event size, maximum gap size across which to merge events, and event types to merge. The gap consideration passes if either the events are close enough together, as specified by the maximum gap size parameter for each of the three individual scores, or the two events span the centromere, with no probes in between them. Adjacent events are smoothed only if they are of identical types, or if they map to the same canonical event according to the characteristics of each scar.

SELECTION

The resulting breakpoints (LST) and calls (TAI and LOH) that comply with each scar's specifications are selected according to each scar's criteria. Scar selection is dependent on whether an event touches the telomere, centromere, or the following event. An event is considered touching the telomere if it is overlapping the region listed in the telomere.txt for its respective chromosome or if it contains the last probe on the arm. An event is considered touching the centromere if it is overlapping the region listed in the centromere.txt file for its respective chromosome or if it contains the last probe on the arm. An event touches the next event if there are no probes in the gap between the two events.

The applied definition of each scar is:

- **Loss of Heterozygosity (LOH)** - number of regions representing one parental allele resulting from a copy number neutral, or a loss, event that is longer than a specified minimum LOH event size, but shorter than the whole chromosome.
- **Telomeric Allelic Imbalance (TAI)** - number of regions with CNV or allelic imbalance longer than the specified minimum TAI event that extends to one of the telomeres but does not cross the centromere.
- **Large-Scale State Transitions (LST)** - number of chromosomal break points between adjacent regions of change in copy number or allelic content longer than a specified minimum LST event size. Adjacent events with a gap less than the maximum LST gap size are merged. State changes at centromeres and telomeres are excluded.

The characteristic event size and gap size for each genomic scar is configurable. A config HRD Parameters file is retained as a TXT file within the VIA server (../VIA Server/Storage/Resources) that can be modified to adjust the

default parameters and refine the scarring performance accordingly. The specific parameters used in calculation of the genomic scars are the minimum event size and the maximum gap size for all three scar types.

HRD Genomic Scar Performance

HRD Genomic Scar Processing was performed on 529 ovarian cancer samples from the Nexus Copy Number TCGA Premier dataset processed in VIA, which had been previously curated to correct for over-segmentation and incorrect ploidy. The combined tally of the genomic scars was compared for 497 samples that had previously reported analyses in [Takaya et al 2020](#) (Table 6). Genomic scar calculations were compared for 96 samples that had either a mutation and/or methylation in either *BRCA1* and/or *BRCA2* and 191 samples that had neither a mutation or methylation in *BRCA1* and *BRCA2* (Table 7, Table 8). As described by Takaya et. al. 2020, samples were defined as being HR deficient if the tally of genomic scars was greater than 63. The status for Nexus Copy Number TCGA Premier samples processed in VIA was defined as being deficient if the tally of genomic scars was greater than 42. Default HRD parameters were used to calculate genomic scar tallies, except TAI gap size, which was set to 0 MB (Takaya, H., Nakai, H., Takamatsu, S. et al. Homologous recombination deficiency status-based classification of high-grade serous ovarian carcinoma. Sci Rep 10, 2757 (2020)).

Table 6. Takaya et al. vs VIA

All Samples (Takaya et al. vs VIA)	
Concordant	414 (83.30%)
Discordant	83 (16.70%)

Table 7. A Mutation/Methylation in at least either *BRCA1/2*

A Mutation/Methylation in at least either <i>BRCA1/2</i>		
HRD Status	Takaya et al. 2020	VIA
Deficient	74 (77.08%)	84 (87.50%)
Proficient	22 (22.92%)	12 (12.50%)

Table 8. No Mutation/Methylation in *BRCA1/2*

No Mutation/Methylation in <i>BRCA1/2</i>		
HRD Status	Takaya et al. 2020	VIA
Deficient	54 (28.27%)	72 (37.70%)
Proficient	137 (71.73%)	119 (62.30%)

The result of the analysis demonstrated high concordance for the automated assessment of genomic scars associated with HRD in VIA as a robust means to determine HR status from technology types producing copy number and B-allele frequency data for application in clinical research oncology, shown in **Figure 13**.

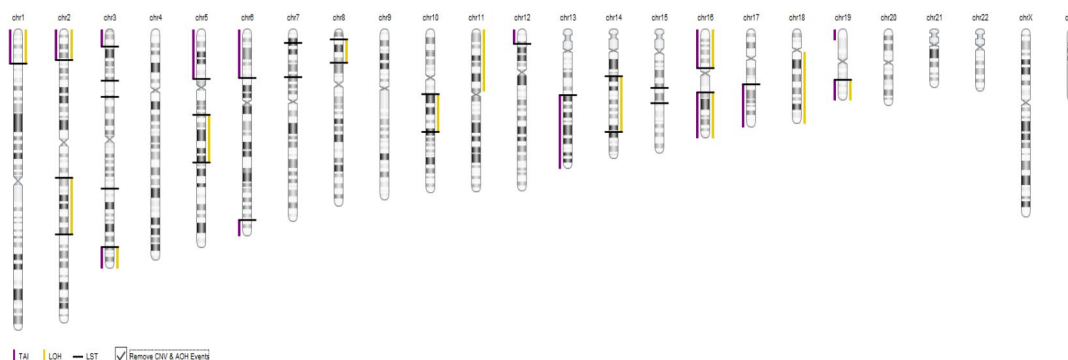


Figure 13. Example visual of HRD genomic scar analysis in VIA

Aneusomy Detection and Analysis

Accurate and sensitive identification of aneuploidy is critical for both the comprehensive characterization of cancers as well as its impact on genetic disease. VIA software offers an innovative method to attain unparalleled sensitivity to detect aneuploidy even at low cell fractions for inclusion in the analysis and reporting workflow. For samples processed with the aneusomy analysis enabled, a table is provided that displays a measurement of estimated copy number for each chromosome (and chromosome arm if selected), an indication of the copy number state designation (Gain/Loss or no call), and a calculation of confidence in the designation.

Aneusomy Calling

Aneusomy calling can be performed on whole chromosomes only or on both whole chromosomes and chromosome arms, which is a configurable parameter with the Processing Settings for the Sample Type. For an aneusomy call to be made, the copy number estimate for the region must exceed the user-modifiable target threshold, specified in terms of the minimum desired aberrant cell fraction (ACF) to detect. Copy number estimates are computed from the median copy number probe value in the region. Calls are made as follows:

$$aneusomy = \begin{cases} LOSS, & m \leq c - a. \\ NORMAL, & c - a < m < c + a \\ GAIN, & m \geq c + a. \end{cases}$$

where m is the median copy number estimate in the region, a is the user specified ACF detection threshold, and c is the chromosome ploidy.

Due to sparseness and high variability in copy number data on the satellited short arm of acrocentric chromosomes, aneusomy calling is not performed on the p arms (indicated as 'NA') and they are excluded from the whole chromosome calculations for acrocentric chromosomes.

Aneusomy Confidence Scoring

A confidence score is included with each aneusomy call ranging from 0 to 1, representing certainty in the existence (or non-existence) of the aneusomy. The confidence value is computed by measuring how significantly the median copy number estimate in the region deviates from the baseline. In the case of an aneusomy call, the baseline is the chromosome ploidy. In the case of a normal call, the baseline lies at the nearest ACF detection threshold. Therefore, when a call is made, a high score indicates high confidence that an aneusomy is present at any level. When a call is not made, a high score indicates high confidence that an aneusomy is not present at the detection threshold.

Deviation is measured by assuming that chromosome (or arm) copy number medians follow a normal distribution. Background medians are collected (excluding the region in question, sex chromosomes, and any other regions that have been called aneusomy) and measure their standard deviation. Subsequently, the p-value is computed assuming a normal distribution parameterized by the baseline as the mean and the measured standard deviation, representing the likelihood of observing the given copy number estimate, in light of the null hypothesis. The final confidence score is the complement of this p-value.

Resolving Conflicts Between Whole Chromosomes and Chromosome Arms

When the user selects the whole chromosome and chromosome arm calling option, it is theoretically possible to observe inconsistencies between whole chromosome and arm calls given that these determinations are made independently. To avoid the influence of chromosome arm level events on the measurement of whole chromosome copy number estimate, a conflict resolving step abstains from making a whole chromosome call in cases where the chromosome has a high disparity between chromosome arm copy number estimates and these estimates are conflicting. High disparity is determined when the copy number estimates for the two chromosome arms fall on opposite sides of the whole chromosome detection threshold, and the two estimates are at least 6 standard deviations apart (as determined by the distribution outlined in the “Confidence Scoring” section) wherein the software refrains from making the whole chromosome call, and instead displays a hashed value.

PERFORMANCE OF ANEUSOMY CALLER WITH OGM DATA (300X EFFECTIVE COVERAGE)

The performance of aneusomy calling was assessed using a set of sample data with simulated events as well as real samples with known aberrations with a 7% ACF aneusomy calling threshold. Whole chromosome and arm sensitivities were 95.5% and 100% in simulations and 100% and 100% in real data. Specificity was >99% and >99% for all simulated VAFs and 98% and 96% in real data.

ANEUSOMY CALLING RESULTS FOR SIMULATED SAMPLE

Datasets were generated that simulated whole chromosome monosomy and trisomy for each autosome at 5%, 10%, and 50% VAF (10%, 20%, 100% ACF) totaling 44 datasets at 300x effective coverage. In addition, eight datasets were generated with chromosome arm level events that are relevant to heme disease analyses (1q del, 1q dup, 7q del, 9q del, 11q del, 17p del, 17q dup, and 20q del) and simulated to the same allele fractions and coverage.

	0.05 VAF		0.1 VAF		0.5 VAF	
	sensitivity	specificity	sensitivity	specificity	sensitivity	specificity
Whole Chromosome (n=44)	0.955	0.999	1	0.995	1	0.999

Chromosome Arm (n=8)	1	1	1	>0.999	NA	0.991
----------------------	---	---	---	--------	----	-------

ANEUSOMY CALLING RESULTS FOR REAL SAMPLES WITH KNOWN ABERRATIONS

A set of samples with whole chromosome and chromosome arm aberrations was used to assess the performance of the aneusomy caller on real samples analyzed with OGM. Expected calls for this sample-set (**Table 9** and **Table 10**) was established from orthogonal assays, and manual curation of the OGM data when genome-wide CNV analysis was not part of the orthogonal results. Aneusomies present in the cell population lower than 5% VAF were excluded from the performance calculations.

Table 9. Aggregated performance of expected calls for real sample dataset

Expected Aneusomy Calls		Sensitivity	Specificity
Total samples n=14	Whole Chromosome (n=8)	100% (8/8)	98.4% (254/258)
	Chromosome arm (N=11)	100% (11/11)	97.1% (437/450)

Table 10. Detailed table of expected results of real sample data

Sample	Sex	Expected Calls	Aneusomy Calls	Comments
221012RR1-01	male	-7	-7, -7p, -7q	
221106ND1	male	+1q, -13	+1q, -13, -13q, -Yp	
B21-2286 (B211039)	male	-2p, +5p, +8, -17p	-2p, +5p, +8, +8p, +8q, -17p, -Y, -Yp, -Yq	
BNGOHH-0000597	female	none	none	
G169_BV_6	female	+1q, +7, -15, -16q, +18q	+1q, +7, +7p, +7q, -8p, +8q, +11p, -15, -15q, -16q, -17p, +18q, +X, +Xp, +Xq	-8p, +8q, +11p, -13q, +13q, -17p large CNVs/many gains on X
HS_CLL_12112_R	female	+3q, +7p, +18	+3q, +7p, -13, -13q, +18, +18p, +18q, -Y, -Yp, -Yq	+7q, +8q, -13q CNVs
LL190091	female	none	+4q	+4q large CNV
OGM24_G2	female	-7q, -8p	-7q, -8p	
OGM33_G2	male	none	none	+4 (2%VAF)
R21-2832	male	none	none	
R23-852	male	none	none	

Observed whole chromosome and chromosome arm aneusomy events called by VIA 7.1 for good quality OGM data of samples with known findings present at >5%VAF

Expected Workflow for Integrating Aneusomy Calls into the Results Table

This section provides guidance on how the aneusomy calls made by the software algorithm can be included in the results table for inclusion in a report.

1. Review the aneusomy calls in the Aneusomy table, which includes details for each chromosome (and chromosome arm if selected) with an indication of an aneusomy call (**GAIN**, **LOSS**, or no call), copy number estimate, and confidence calculation, as represented in the partial table example below.

	chr1	chr2	chr3	chr4	chr5	chr6	chr7	chr8
Aneusomy Call	-	-	-	-	-	-	-	GAIN
Copy number estimate	1.99	1.98	1.98	1.98	1.99	1.99	2.01	2.17
Confidence	> 0.99	> 0.99	> 0.99	> 0.99	> 0.99	> 0.99	> 0.99	> 0.99

2. Navigate to the desired chromosome by selecting the column header.
3. Manually create whole chromosome or arm calls from the **Probes** track/view in the **Tracks** tab by:
 - a. Whole Chromosome call- Right clicking on the **Probes** track and selecting **Whole Chromosome Gain** or **Whole Chromosome Loss** will set the ISCN for an aneusomy appropriately, seen in

Figure 14.

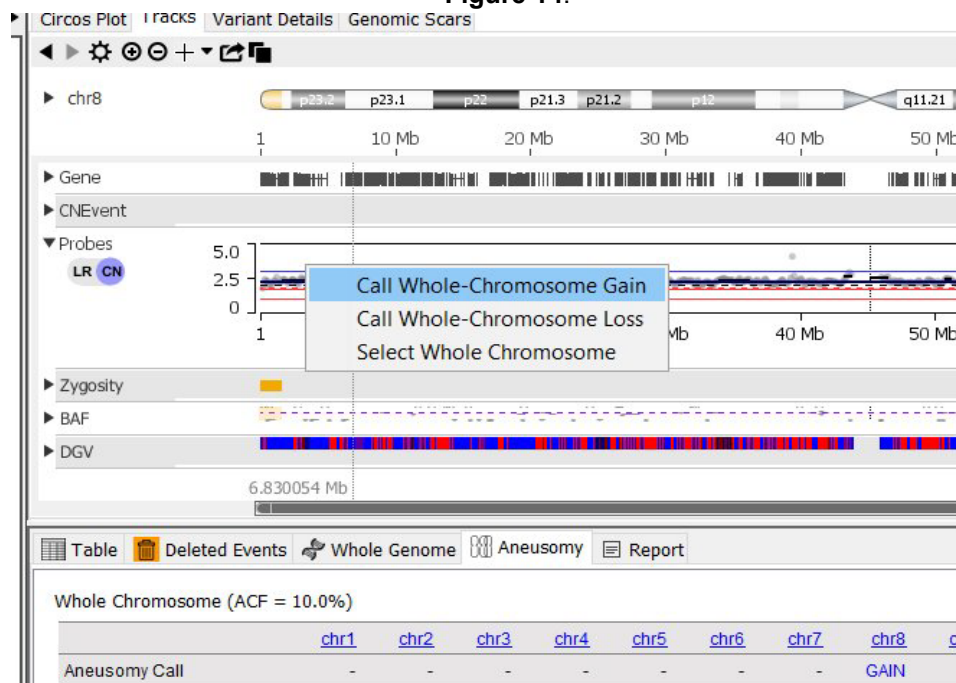


Figure 14. Whole Chromosome call

- b. Arm call – Use the selection tool to select the arm region (**Figure 15**); right click to make the call.

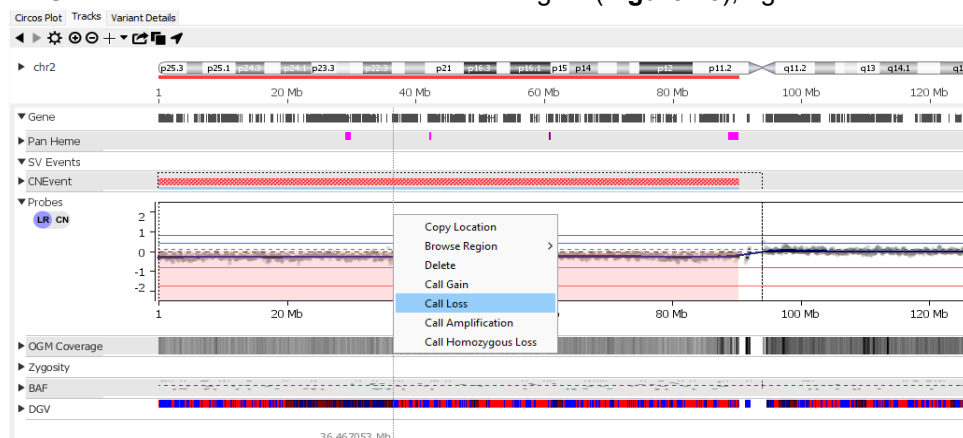


Figure 15. Arm call

- This manual operation will rectify fragmented CNV calls to generate a single call visible in the **Table** and **Tracks** views. Subsequently, the aneusomy calls will be included in all visuals and the **Event Table** for reporting and will be included within the CN events database.
- Upon event creation, a prompt will display (when available) inquiring whether to run auto-classification of the aneusomy event to classify according to the defined decision tree.

Guidance for Interpreting Aneusomy Confidence Scores

CONSIDERATION OF SAMPLE QUALITY AND GENOME COMPLEXITY

The quality of the sample should be considered when reviewing the table of aneusomy calculations (see **Figure 16**). Lower quality samples will display a high variance of CN Estimate values for each chromosome from 2.0 (disomy) and caution should be used to extract a higher sensitivity to detect aneusomies. In VIA, good quality samples will typically have a low CN Quality value indicating more reliable data. The CN Quality values are only used to compare relative values between samples within a platform. For OGM data, good quality samples will have a CN Quality value < 0.01 in VIA, while additional metrics to assess sample quality are provided in *Bionano Solve Theory of Operation Structural Variant Calling* (CG-30110).

	chr1	chr2	chr3	chr4	chr5	chr6	chr7	chr8	chr9	chr10	chr11	chr12	chr13	chr14	chr15	chr16	chr17	chr18	chr19	chr20	chr21	chr22
Aneusomy Call	-	-	-	-	-	-	-	-	-	-	-	-	-	-	-	-	-	-	LOSS	-	-	LOSS
Copy number estimate	2.02	2.00	2.00	1.99	2.01	2.01	1.97	1.99	1.99	1.99	2.00	2.03	2.00	2.04	2.03	1.98	2.03	1.95	1.91	1.96	1.95	1.90
Confidence	0.96	> 0.99	> 0.99	0.98	0.98	0.98	0.92	0.99	0.99	0.98	> 0.99	0.94	> 0.99	0.90	0.93	0.96	0.94	0.77	> 0.99	0.90	0.82	> 0.99

Figure 16. Example of Aneusomy results for a sample with failing QC metrics

- It is common for a sample with polyploidy or gross genomic complexity to not center on the baseline properly, resulting in high CN Quality and BAF Quality values and spurious CN Estimate values in the Aneusomy table. Recentering the genome on a region of disomy will reset and correct the Aneusomy CN Estimate values and calls accordingly.

GUIDANCE TO IDENTIFY POSSIBLY TRUE UNCALLED EVENTS AT VERY LOW VAF (ACF)

Aneusomies with a VAF (ACF) lower than the specified threshold are not expected to be called. However, it is possible to detect and evaluate very low mosaic true positives using the following approaches as guidance:

- Disomic chromosome calls that are approaching the specified aneusomy calling threshold will have low confidence values and will be highlighted in the aneusomy table (<95%):


Q-arm (ACF = 10.0%)

	chr1
Aneusomy Call	-
Copy number estimate	2.08
Confidence	0.82

- Perform manual inspection of CN Estimate values in the aneusomy table to identify those distinguishable from other chromosomes. The effectiveness of this subjective review will vary by analyst.

	chr15	chr16	chr17	chr18	chr19	chr20	chr21	chr22
Aneusomy Call	-	-	-	-	-	-	-	-
Copy number estimate	2.00	1.99	1.99	2.06	1.99	1.98	1.99	1.99
Confidence	> 0.99	> 0.99	> 0.99	0.92	> 0.99	> 0.99	> 0.99	> 0.99

- Using the aneusomy table as an aid, perform manual inspection of the whole genome Copy Number and BAF plot to identify observed shifts in the moving average implying the presence of a low level aneusomy (see example below). The aneusomy table data and whole genome plots combined assist to inform legitimacy of suspected very low VAF aneusomy events but the effectiveness of this subjective review will vary by analyst.

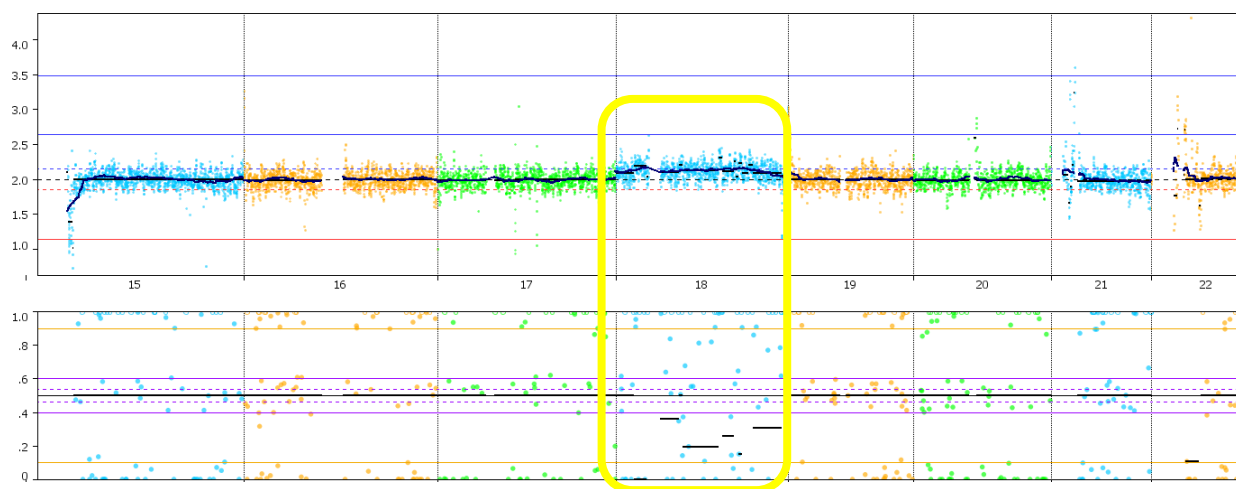

- Very low VAF aneusomy can be added to the results table using the same workflow for integrating aneusomy calls described above.

FACTORS TO CONSIDER THAT COULD LEAD TO FALSE POSITIVE ANEUSOMY CALLS

The aneusomy caller, which seeks to call numerical aberrations affecting the whole of the chromosome or chromosome arm, is susceptible to false positive calls due to the presence of large segmental copy number changes that are biological in origin. False positive calls are possible for both whole chromosomes and arm calls due to the following situations:

- Whole chromosome aneusomy calls can be impacted by multiple large segmental CNVs across both chromosome arms.
- Chromosome arm aneusomy calls can be impacted by large segmental CNVs on the chromosome arm.

Sample quality affects false positives as well. Lower quality samples will display a high variance of CN Estimate values from 2.0, with potential for high percentage of chromosomes to display low confidence values.

Automatic Pre-Classification Decision Trees

To specify rules for automatic pre-classification, one must create the rules using a specified syntax. Failure to follow the syntax carefully will result in errors and the automated classification may not work. Care also must be taken to make sure all parentheses and curly braces match. The functions used for the decision tree rules are case-sensitive so attention must be given to this as well. Specific keywords and syntax used to create the decision tree rules are detailed in this section. **NOTE:** Case, spacing, and use of quotes are important so please make note of that in the examples. Any quotes must be ASCII quotes. Type in quotation marks directly into the query field rather than copying and pasting from other documents as often these use smart quotes rather than ASCII quotes. If quotes are not ASCII quotes, the query will not work.

Rules typically follow the IF, THEN, ELSE statement syntax with use of logical operators **AND** and **OR**. For example:

```
IF {A > 10} THEN {X = "yes"}
IF {(A < 10) AND (B > 10)} THEN {X = "no"}
```

There are also other functions used to test various attributes of the event in question.

An example of a rule:

To classify as **Likely Benign**, any copy number loss event that does not overlap any gene in the **OMIM Morbid** track, the following rule can be specified:

```
IF {CN_EVENT_KIND IS CN_LOSS OR CN_EVENT_KIND IS HOMOZYGOUS_COPY_LOSS}
THEN {
  IF {!(OVERLAP(OMIM Morbid) > 0 )}
  THEN {CLASSIFY("Likely Benign")}
}
```

Or it can be written more concisely using CASE statements:

```
CASE {CN_EVENT_KIND IS CN_LOSS OR CN_EVENT_KIND IS HOMOZYGOUS_COPY_LOSS}
{ CASE{!(OVERLAP(OMIM Morbid) > 0 )} {CLASSIFY("Likely Benign")} }
```

NOTE: there are examples of decision tree scripts in the last section of this document.

Keywords and Syntax For the Decision Tree Rules

TEXT STYLE

All functions must be in uppercase, or the function will not work.

```
IF {A > 10} THEN {X = "yes"}  works
if {A > 10} Then {X = "yes"}  will not work
```

COMMENTING IN THE SCRIPT

Two forward slashes (//) can be used to add comments to the scripts and text but this will not be interpreted as part of the code. For example:

```
CASE {CN_EVENT_KIND IS CN_LOSS OR CN_EVENT_KIND IS HOMOZYGOUS_COPY_LOSS}
{ CASE{!(OVERLAP(OMIM Morbid Phenotypes) > 0 )} {CLASSIFY("Likely Benign")} } // calls copy losses that do not
overlap with the OMIM Morbid Phenotypes track as "Likely Benign"
```

AUTOCOMPLETE

So that users do not have to remember all various functions and possible values, the decision tree editor has a predictive text/autocomplete feature where as one starts typing, possible values are displayed. Use the mouse to select a value by double clicking on it, using the **TAB** key to select the first item or use the arrow keys to highlight the term and click **Enter**. **NOTE:** Autocomplete only works if the user is typing using ALL CAPS.

Typing an uppercase **C** brings up the list seen in **Figure 17**. If a lowercase **c** is entered, no options are displayed for completing the statement. **Figure 18** is an example of looking for CN events (the **Ctrl-Space** command can also be keyed in to see a list of all possible keywords).

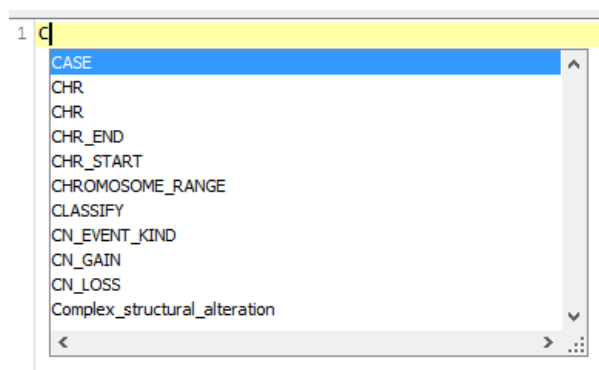


Figure 17. Autocomplete



Figure 18. Looking for CN Events

OPERATOR BASICS

When the value being evaluated is numerical, the following operators are used:

- > Greater than
- < Less than
- >= Greater than or equal to

`!=` Not equal to

`<=` Less than or equal to

`==` Equal to (**NOTE:** must use two equal signs in succession)

When the value being evaluated is a text string, the following operators are used:

`==` Equals

`!=` Does not equal

Basic Functions

VAL (USING VARIABLES)

It is possible to define and set a variable at the beginning of the script to be used throughout. This makes it easy to change various thresholds that may be defined by changing the value of the variable rather than searching for and changing the constant in many different locations.

The keyword **val** is used to define a variable, e.g.,

```
val percent_overlap = 0.50
```

```
val similarity_threshold = 0.85
```

```
val listSampleTypes = List("CytoScan HD," "Affymetrix OncoScan")
```

Then within the script the variables can be used instead of the actual values:

```
CASE {OVERLAP("OMIM Morbid Phenotypes") > percent_overlap} {CLASSIFY("Likely Pathogenic")}
CASE {SIMILARITY("DECIPHER Syndromes") > similarity_threshold} {CLASSIFY("Pathogenic")}
CASE {SIMILARITY ("DECIPHER Syndromes Gains")> similarity_threshold } {CLASSIFY("Pathogenic")}      // Classify
pathogenic based on similarity to cases in DECIPHER Syndromes
CASE {SIMILARITY ("ClinGen Dosage Sensitive Map Triplosensitivity Pathologic Regions")> similarity_threshold }
{CLASSIFY("Pathogenic")}
CASE { (SIMILAR_CASES("ClinGen Postnatal Gains Uncertain Significance", similarity_threshold) >= 3) OR (SIMILAR_CASES
("ClinGen Postnatal Gains Pathogenic", similarity_threshold) >=3 ) OR (SIMILAR_CASES ("ClinGen Postnatal Gains
Likely Pathogenic", similarity_threshold) >=3 )} {CLASSIFY("VUS")}
```

If the numbers 0.5 and 0.85 are used in many places in the script and adjustments are desired, change the numbers in one spot. To use 0.89 instead of 0.85 as the similarity threshold, change the value of the variable once at the beginning of the script.

IF {} THEN {} ELSE {}

If the condition being evaluated is true do one thing and if it is false do something else.

```
IF { condition is true} THEN {do this} ELSE {do that}
IF {OVERLAP("OMIM Morbid Phenotypes") > 0} THEN {CLASSIFY("Possible reportable event")} ELSE {CLASSIFY("VUS")}
```

If the event overlaps with a region in the **OMIM Morbid Phenotypes** track, then classify the event as “Possible reportable event” or else classify it as “VUS.”

The following uses of IF/THEN statements are NOT supported:

- Combining IF/THEN statements in a decision tree script without ELSE clause:

```
IF {ANY_SNP_EVENT_KIND} THEN { CLASSIFY("VUS") }

IF {EVENT_SIZE > (1 Mb)} THEN { CLASSIFY("Pathogenic") }
```

- Nest IF statements:

```
IF {ANY_SNP_EVENT_KIND} THEN {CLASSIFY("VUS")} ELSE {IF {EVENT_SIZE > (1 Mb)} THEN {CLASSIFY("Pathogenic")}}
```

AND/OR

AND or **OR** can be used to check that either multiple conditions are true or that one of multiple conditions is true. For example:

```
IF {(grade > 1) AND (grade < 6)} THEN {school = "elementary"} ELSE {school = "not elementary"}

IF {(grade == 6) OR (grade == 7) OR (grade == 8)} THEN {school = "middle"} ELSE {school = "not middle"}
```

Example usage:

```
IF { (OVERLAP("OMIM Morbid Phenotypes") < 0) AND (OVERLAP("RefSeq Exons")==0) } THEN {CLASSIFY("Likely Benign")}
ELSE {CLASSIFY("Unclassified")}
```

CASE {} {}

This replaces multiple **IF/THEN** statements and is the preferred way to write statements as it is more concise. The first set of parentheses has a condition being tested, if the condition is true, the second set contains the action. Additional CASE statements can be added to the line. Conditions are evaluated in sequence until a stopping/terminating action is encountered (e.g., CLASSIFY). Additional CASE statements can be nested within the second set of curly braces (the action item). If a true condition is not met, evaluations move on to the next rule.

```
CASE {if grade equals "kindergarten"} {then in elementary school}
CASE {if grade <= 5} {then in Elementary school}
CASE {if grade <= 8} {then in Middle School}
CASE {if grade <= 12} {then in High School}
```

If in the example above, the grade happens to be 13, then none of the cases will be true and the next rule will be evaluated. If the grade happens to be 7, then the first three cases will be evaluated and the 3rd case will result in *true*. The next case will not be evaluated, and the next rule will be evaluated.

Example: if the DGV score is greater than .98 then classify as “Benign.”

```
CASE {SCORE(DGV) > .98} {CLASSIFY("Benign") }
```

If the DGV score is greater than .98 then classify as Benign. If not, then evaluate the second case to see if the score is greater than .85 and, if so, then classify as “Likely Benign.”

```
CASE {SCORE(DGV) > .98} {CLASSIFY("Benign")} }
CASE {SCORE(DGV) > .85} {CLASSIFY("Likely Benign")} }
```

NOTE: In the nested CASE statement below, if the event size is less than 100 then the next CASE is evaluated. If the event is a copy number event, then it is classified as “Likely Benign;” if it is not a copy number event or not less than 100, the next rule would be evaluated.

```
CASE {EVENT_SIZE <100} { CASE {CN_EVENT_KIND IS COPY_GAIN} {CLASSIFY("Likely Benign")}}
```

Classifying an Event

CLASSIFY()

Classifying an Event allows for termination of the script. The value that goes inside the parentheses must be a text string in quotes and it must match one of the pre-defined classification values set up through the VIA Administrator. For example, the VIA Administrator has defined the following classification values: “Pathogenic” “Likely Pathogenic,” “Benign,” “Likely Benign,” “VUS.”

Example usage:

```
CASE {OVERLAP("OMIM morbid")} {CLASSIFY("Likely Pathogenic")}
```

NOTE: If the event overlaps with an event in the **OMIM Morbid Phenotypes** track, then the event is probably pathogenic so classify this as “Likely Pathogenic.”

The syntax above is valid (the value passed into “CLASSIFY” is a defined classification term and is enclosed in quotes). The following are examples of invalid syntax usage:

```
CASE {OVERLAP("OMIM morbid")} {CLASSIFY("LIKELY PATHOGENIC")}
```

NOTE: The value passed into CLASSIFY is case-sensitive and must exactly match the pre-defined classification values. Here the value passed in is in all uppercase whereas the defined pre-classification values are not (defined value: Likely Pathogenic).

```
CASE {OVERLAP("OMIM morbid")} {CLASSIFY(Maybe Pathogenic)}
```

NOTE: The value passed in to CLASSIFY must match exactly the pre-defined classification values. “Maybe Pathogenic” is not one of the pre-defined classification values. Also, note that quote marks are required.

```
CASE {OVERLAP("OMIM morbid")} {CLASSIFY(Likely Pathogenic)}
```

NOTE: Value passed into CLASSIFY is a pre-defined classification value, but it is not enclosed in quotes, which is required.

Evaluating the Type of Event (Copy number gain, AOH, SNV, etc.)

These functions assess the type for the event. They evaluate whether the event is a copy number event and what type of copy number event (single copy gain, high copy gain, ...) or whether the event is an allelic event (AOH, total allelic loss, ...).

ANY_CN_EVENT_KIND, ANY_SNP_EVENT_KIND, ANY_SEQVAR_EVENT_KIND

Looks for any copy number (ANY_CN_EVENT_KIND), allelic (ANY_SNP_EVENT_KIND), or sequence variant (ANY_SEQVAR_EVENT_KIND) event.

NOTE: Prior to version 7.0 ANY_SV_EVENT_KIND is used for sequence variants. Version 7.0 onwards use ANY_SEQVAR_EVENT_KIND for sequence variants.

Example usage:

```
CASE {ANY_CN_EVENT_KIND}
{ CASE{SCORE(DGV) > .95} {CLASSIFY("Likely Benign")} }
```

If the event is any copy number event and the DGV score is greater than 0.95, then classify as “Likely Benign.”

CN_EVENT_KIND IS, SNP_EVENT_KIND IS, SEQVAR_EVENT_KIND IS

NOTE: Prior to version 7.0 SV_EVENT_KIND IS was used for sequence variants. Version 7.0 onwards use SEQVAR_EVENT_KIND IS for sequence variants.

Checks for the type of copy number event (e.g., high copy gain), SNP event (e.g., AOH), or sequence variant event (e.g., SNV, Insertion, Deletion).

Example usage:

```
CASE {CN_EVENT_KIND IS HIGH_COPY_GAIN} {...}
CASE {SV_EVENT_KIND IS STRUCTURAL_ALTERATION} {...}
```

Valid values for **CN_EVENT_KIND**:

HIGH_COPY_GAIN
CN_GAIN
CN_LOSS
HOMOZYGOUS_COPY_LOSS

Valid values for **SNP_EVENT_KIND**:

AOH
ALLELIC_IMBALANCE
HETERO_UPD_FATHER
HETERO_UPD_MOTHER
ISO_UPD_FATHER
ISO_UPD_MOTHER

Valid values for **SEQVAR_EVENT_KIND**:

SNV, Indel, ...

There are numerous possible values for SEQVAR_EVENT_KIND. To see the complete list of possible values, go to the **Admin->Sample Types->Sample Review Preferences->Filter** tab. Any values listed in the “Seq. Var. Events” “Show Events of Type” section are permitted, as shown in **Figure 19**.

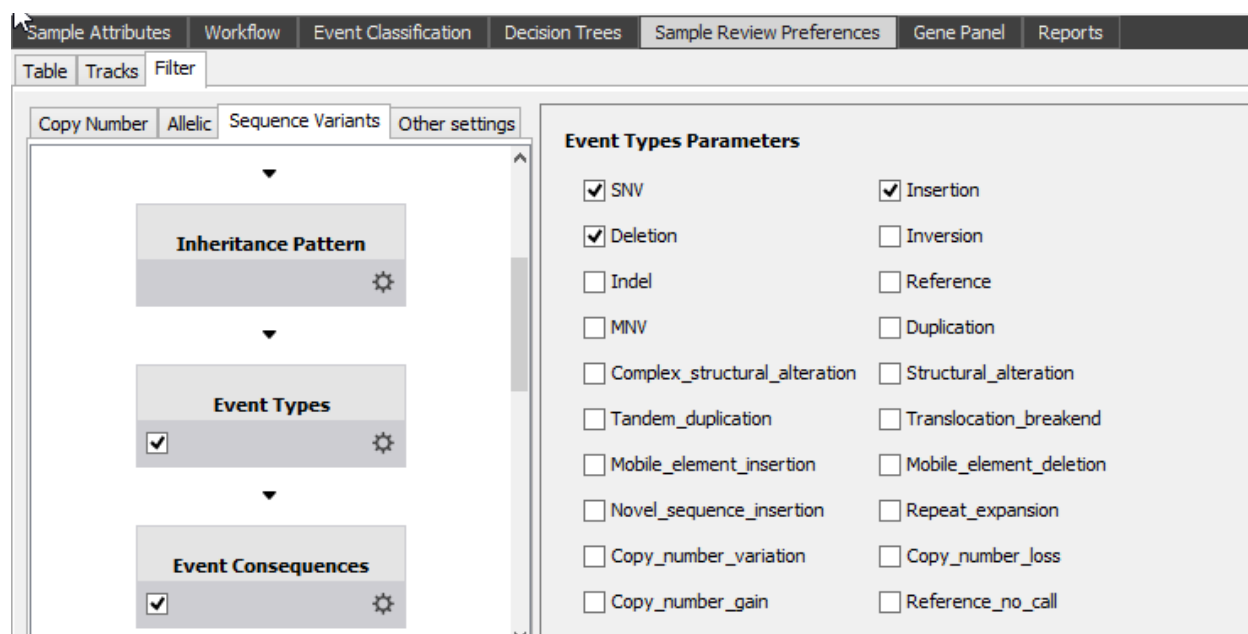


Figure 19. Screenshot of partial list of values in the **Event Types** filter.

The negation can be used with this function to look for any event other than a single event type. For example, if one wants to check to see if a CN event is a high copy gain and nothing else, one can use the negation (!):

```
CASE {(ANY_CN_EVENT_KIND) AND !(CN_EVENT_KIND IS HIGH_COPY_GAIN)} {...}
```

If the event is a copy number event other than high copy gain, then...

NOTE: In the statement above, users must also first check to see if the event is a copy number/allelic/sequence variant event using the functions beginning with “ANY_.” For example, for copy number, only checking to see that it is not a high copy gain will return any other event (including allelic events) as true.

Evaluating Event Size and Location

EVENT_SIZE

Specifies the length of an event (gain, loss, etc.). Units can be in bp, Kb, Mb, and Gb. The following all equal 1,000,000,000 bp:

1000000000

1000000.Kb

1000.Mb

1.Gb

NOTE: To use Gb, Mb, Kb, a period must be placed between the number and the unit. Specifying as 1000 Mb is incorrect. It must be specified as 1000.Mb. When specifying in bp, no unit is needed.

Allowed operators:

- > Greater than
- < Less than
- >= Greater than or equal to
- != Not equal to
- <= Less than or equal to
- == Equal to (**NOTE:** must use two equal signs in succession)

Example usage:

```
CASE {EVENT_SIZE >= 5000} {...}
CASE {EVENT_SIZE >= 5.Kb} {...}
```

CHR IS

Determines on which chromosome the event is located. The function returns a “true” or “false.”

Example usage:

```
CASE {CHR IS Chr4} {do this}
```

Syntax for specifying chromosomes: use uppercase for first letter or all lowercase with chromosome numbering ranging from 1 to 24 and also use X and Y, if desired.

For example, the following are all valid:

```
Chr3, chr3, ChrX, chrX, ChrY, chrY, Chr23, chr23
```

POSITION()

Refers to the location of an event and returns a “true” or “false.” Arguments (locations on one or more chromosomes, either a base pair location range or a chromosomal location range) are passed in to the POSITION() function separated by commas. An optional minimum overlap value can be passed in as the last item in the argument list and if so, the function will then test to see if at least the specified percentage of the event is overlapping with one of the specified regions. If no minimum overlap value is explicitly specified, it is assumed to be 1.0 (i.e., 100% of the event must lie within one of the specified ranges). The minimum overlap argument must be a number between 0 and 1.0.

If the minimum overlap value is not specified, an unlimited number of chromosome ranges can be passed in. If minimum overlap is specified, up to 12 chromosome ranges can be passed in. **NOTE:** The ranges specified are all inclusive (i.e., if the range specified is 10000->20000, then the function will return true if the event is from 10000 to 11000 and if the event is from 17000 to 20000).

Example usage without specifying a minimum overlap:

```
CASE {POSITION(10000 -> 20000, Chr1 -> Chr3, Chr5::CHR_START -> 500000)} {do this }
```

NOTE: If the entire event (100%) falls in any one of the following ranges, then it is true:

- between and including 10,000 bp and 20,000 bp on any chromosome or
- anywhere on chromosomes 1, 2, or 3, or
- on chromosome 5 between and including the start of the chromosome and 500,000 bp

Example usage specifying a minimum overlap:

```
CASE { POSITION (CHR4::CHR_START -> 10.Mb, Chr10::10000000 -> 15000000, 0.6) } { do this }
```

NOTE: If at least 60% of the event falls in any one of the specified ranges, then the function will evaluate to true.

NOTE: For events that are much larger than the region, even if the specified region is completely covered/overlapping with the event, the function may not return as true since the specified percentage of the event must overlap the region. For example, if the region is 200 bp and the event is 600 bp and the event completely covers the region, the following function will not return as true since 60% of the event (360 bp) cannot overlap the region as the region is too small (200 bp).

```
CASE { POSITION (Chr10::10000000 -> 10000200, 0.6) } { do this }
```

Location units can be in bp, Kb, Mb, and Gb. The following all equal 1,000,000,000 bp:

- 1000000000
- 1000000.Kb
- 1000.Mb
- 1.Gb

NOTE: To use Gb, Mb, Kb, a dot must be placed between the number and the unit. Specifying as 1000 Mb is incorrect. It must be specified as 1000.Mb. When specifying in bp, no unit is needed. Example:

```
CASE { POSITION(10.Kb -> 2.Mb, Chr1 -> Chr3, Chr5::CHR_START -> 500.Kb, Chr10::CHR_START->200000) } {do this}
```

CHR_START, CHR_END

Refers to the start and end positions of a chromosome. **CHR_START** refers to the bp location of the start of the chromosome and **CHR_END** refers to the bp location of the end of the chromosome. These can only be used in conjunction with the **POSITION** function.

Example usage:

```
CASE { POSITION(10000 -> CHR_END) } {do this }
```

NOTE: If the event falls within 10,000bp and the end of the chromosome on any chromosome, then:

```
CASE { POSITION(CHR_START -> 200000) } {do this}
```

NOTE: If the event falls within the start of the chromosome and 200,000bp on any chromosome, then:

```
CASE { POSITION(CHR_START -> 200000, 300000 -> CHR_END) } {do this}
```

NOTE: If the event falls within the start of the chromosome and 200,000bp, or it falls within 300,000bp and the end of the chromosome, then:

```
CASE { POSITION(Chr3::100000 -> CHR_END) } { ... }
```

NOTE: If the event falls on chromosome 3 within 100,000bp and the end of the chromosome, then:

```
CASE { POSITION(Chr3::CHR_START -> 200000) } { ... }
```

NOTE: If the event is on chromosome 3 and falls within the start of the chromosome and 200,000bp, then:

```
CASE { POSITION(Chr3::CHR_START -> 200, Chr4::300 -> CHR_END) } { ... }
```

NOTE: If the event is on chromosome 3 and falls within the start of chromosome 3 and 200,000bp or if the event is on chromosome 4 and falls within 300,000bp and the end of chromosome 4, then:

Comparing Events to Region lists and Evaluating Similarity Scores

EVENT_OVERLAP()

Replaces the function OVERLAP().

Tests to see how much of the event overlaps with regions in a specified region list and how the overlap percentage compares to a specified threshold. The input to EVENT_OVERLAP() is a text string that is the name of an annotation list/track. The EVENT_OVERLAP function will then determine if the event in question overlaps with (is covered by) one or more regions in the list and will return a decimal number indicating the percent of the event overlapping the region. For example, if the event overlaps with a region in the specified region list by 31% (i.e., 31% of the event overlaps with the specified region), the EVENT_OVERLAP function will return a value of “.31”. Then the percent overlap can be compared to a given percent using operators to return either TRUE or FALSE. See examples in **Figure 20, Figure 21, and Figure 22.**



Figure 20. Example 1: EVENT_OVERLAP (OMIM Morbid Phenotypes) will evaluate to 1 (100%)



Figure 21. Example 2: EVENT_OVERLAP (OMIM Morbid Phenotypes) will evaluate to 0.8 (80%)

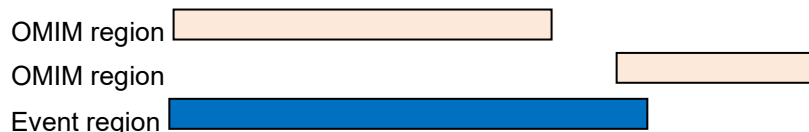


Figure 22. Example 3: EVENT_OVERLAP (OMIM Morbid Phenotypes) will evaluate to 0.9 (90%) since 90% of the event is covered by regions overlapping the event.

When multiple regions are overlapping the event, the coverage will be added or summed. It is possible for this to be greater than 1 if many regions stack up over the sample event.

Allowed operators:

- > Greater than
- < Less than
- >= Greater than or equal to
- != Not equal to
- <= Less than or equal to
- == Equal to (**NOTE:** Must use two equal signs in succession)

A new optional, padding/flanking region parameter has been introduced in version 7.0. This parameter is used to specify a base pair padding to the **Event Region** in both directions when calculating the overlap. The padding region can be specified without units to indicate base pairs or with units, Kb, Mb, and Gb.

Example usage:

```
CASE {EVENT_OVERLAP("OMIM Morbid Phenotypes") > .51} {CLASSIFY("Likely Pathogenic")}
CASE {EVENT_OVERLAP("OMIM Morbid Phenotypes", 2000000) > .51} {CLASSIFY("Likely Pathogenic")}
CASE {EVENT_OVERLAP("OMIM Morbid Phenotypes", 2 Mb) > .51} {CLASSIFY("Likely Pathogenic')}
```

If more than 51% of the event overlaps a region in the **OMIM Morbid Phenotypes** track, then classify the event as “Likely Pathogenic.”

REGION_OVERLAP()

Tests to see how much of a region in a specified region list overlaps with the event in question and how the overlap percentage compares to a specified threshold. The input to REGION_OVERLAP() is a text string that is the name of an annotation list/track. The **REGION_OVERLAP** function will then determine if a region overlaps with (is covered by) the event in question and will return a decimal number indicating the percent of the region overlapping the event. For example, if a region in the specified region list overlaps with the event by 31% (i.e., 31% of the region is covered by the event), the **REGION_OVERLAP** function will return a value of “.31.” Then the percent overlap can be compared to a given percent using operators to return either TRUE or FALSE. See examples in **Figure 23**, **Figure 24**, and **Figure 25**.

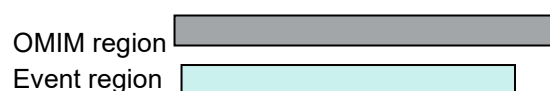


Figure 23. Example 1: REGION_OVERLAP(OMIM Morbid Phenotypes) will evaluate to 0.8 (80%).



Figure 24. Example 2: REGION_OVERLAP(OMIM Morbid Phenotypes) will evaluate to 1 (100%).



Figure 25. Example 3: REGION_OVERLAP (OMIM Morbid Phenotypes) will evaluate to 1 (100%).

In cases where multiple regions are overlapping the event, each region's overlap with the event is calculated separately and then the largest overlap value is returned. Since one region is completely covered by the event (100%) and the other is only covered 20% by the event (value is .2), the larger value, 1 (100%) is returned. This value cannot be greater than 1.

Allowed operators:

- > Greater than
- < Less than
- >= Greater than or equal to
- != Not equal to
- <= Less than or equal to
- == Equal to (**NOTE:** Must use two equal signs in succession)

Example usage:

```
CASE {REGION_OVERLAP(OMIM Morbid Phenotypes) > .80} {CLASSIFY("Likely Pathogenic")}
```

NOTE: If more than 80% of a region in the **OMIM Morbid Phenotypes** track overlaps the event, then classify the event as "Likely Pathogenic."

SIMILARITY()

Tests to see if the event is like regions (of the same class) in a specified region list. The input to SIMILARITY() is a text string that is the name of an annotation list/track. The **SIMILARITY** function will then determine if the event in question is like a region in the list and will return a decimal number indicating the percent similarity. It only looks at those events that are of the same class/type (e.g., if the event in question is a gain, then the function only looks at gain events in the region list). For example, if the event is like a region in the specified region list by 70%, the SIMILARITY function will return a value of "0.7". Then the similarity can be compared to a given percent using operators to return either TRUE or FALSE.

Similarity to another event is based on a ratio of overlap of the common region between the two events and the entire length encompassed by the two events. The percent similarity is defined as a/b where a is the event in question and b is the region covered by both the region being compared as well as the event. **Figure 26** demonstrates this concept.

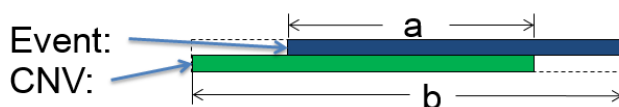


Figure 26. Similarity

Allowed operators:

- > Greater than
- < Less than
- >= Greater than or equal to
- != Not equal to
- <= Less than or equal to
- == Equal to (**NOTE:** Must use two equal signs in succession)

Example usage:

```
CASE {SIMILARITY("DECIPHER Syndromes") > 0.5} {CLASSIFY("Pathogenic")}
```

NOTE: If the event is like an event in the DECIPHER Syndromes list by more than 50%, then classify the event as "Likely Pathogenic".

SIMILAR_CASES()

Tests to see if the event in question is like events in a region file in the system and if the events meet or exceed a similarity threshold also passed in as an argument. The input to SIMILAR_CASES() is a text string that is the name of a region file and a number indicating the similarity threshold. The number of similar events that pass that criteria are then compared to a given number using operators to return either TRUE or FALSE.

SIMILAR_CASES (region file, similarity threshold): This looks at events in a region file like the event in question at or exceeding a similarity threshold also passed in as an argument. The number of similar events that pass that criteria are then compared to a given number using operators to return either TRUE or FALSE. For example, **SIMILAR_CASES**("ClinGen Postnatal Gains Pathogenic", 0.9) looks at events in the "ClinGen Postnatal Gains Pathogenic" region file that are similar to the event in question by at least 90% (see section on *SIMILARITY*() for details on calculation of similarity). Next, the function makes sure that the number of cases meeting the similarity requirement meets the user defined minimum number of cases. E.g.,

```
CASE {(SIMILAR_CASES("ClinGen Postnatal Gains Pathogenic", 0.9)) >= 6 } { ... }
```

NOTE: If the "ClinGen Postnatal Gains Pathogenic" region file contains events that are like the event in question by at least 90% and if there are at least six such events in the region file, then move to the next step.

SCORE()

A comparison of specific properties/functions to the event in question is made to generate a score. It is used in conjunction with other functions to provide the number of previous cases with similar events, the DGV score, and the SAP score.

Keywords/functions that can be passed to the Score() function: **DGV, PREVIOUS_SIMILAR_CASES, SAP**

PREVIOUS_SIMILAR_CASES()

This function looks at previous cases to find those that are like the case under review based on criteria passed into the function. It only works for copy number events. This function must be used in conjunction with the **SCORE** function in the following manner:

SCORE(PREVIOUS_SIMILAR_CASES (classification, similarity threshold, list of Sample Types, earliest processing date, include duplicate samples))

PREVIOUS_SIMILAR_CASES looks at previous cases to return the number that matches the classification that is passed into the function and that meet or exceed a similarity threshold also passed in as an argument. Additional parameters that can be passed in include a list of Sample Types to limit the search to only those Sample Types and the earliest processing date to limit the search to only those samples processed on or since the specified date. Only the classification parameter is required and others are optional, but parameters must be specified in this order: classification, similarity threshold, list of Sample Types, earliest processing date. If a similarity threshold is not specified in the function, a default of 0.9 is used. The “include duplicate samples” argument is an optional boolean field. If it is not specified, then the default behavior is to exclude duplicate samples.

The number of similar events that pass that criteria are then compared to a given number using operators to return either TRUE or FALSE via the **SCORE** function. For example:

SCORE(PREVIOUS_SIMILAR_CASES(“Benign,” 0.95) > 2) looks at past cases in the database that have been classified as “Benign”. It then identifies those that are similar to the event in question by at least 95% (see section on *SIMILARITY*() for details on calculation of similarity). Next, the function makes sure that the number of cases meeting the similarity requirement meets the user defined minimum number of cases.

The list of Sample Types is passed in as a list using the following syntax:

```
List("Affymetrix CytoScan 750K - Postnatal", "Affymetrix CytoScan HD - Postnatal", "Illumina CytoSNP-850K - Postnatal")
```

E.g., **SCORE**(**PREVIOUS_SIMILAR_CASES**("Benign", 0.95, List("Affymetrix CytoScan 750K - Postnatal", "Affymetrix CytoScan HD - Postnatal", "Illumina CytoSNP-850K - Postnatal"))) >= 6 } { **CLASSIFY**("Benign")}

NOTE: A variable can also be specified for the list of Sample Types at the beginning of the decision tree script and then the variable can be used in the script instead. For example:

```
val limitSampleTypes = List("Affymetrix CytoScan 750K - Postnatal", "Affymetrix CytoScan HD - Postnatal", "Illumina CytoSNP-850K - Postnatal")
```

E.g., **SCORE**(**PREVIOUS_SIMILAR_CASES**("Benign", 0.95, limitSampleTypes)) >= 6 } { **CLASSIFY**("Benign")}

NOTE: To also specify which samples to match based on processing date, an earliest processing date can be specified. This means that the function will only look at samples processed on or after the specified date.

E.g., **SCORE**(**PREVIOUS_SIMILAR_CASES**("Benign", 0.95, limitSampleTypes, "2018-01-01")) >= 6, true } { **CLASSIFY**("Benign")}

If the database contains events that have previously been classified as “Benign” that are similar to the event in question by at least 95%, that are of Sample Types specified by the variable “limitSampleTypes”, where samples were processed on or after January 1, 2018, with duplicated samples excluded and if there are at least six such

similar past cases, then classify the current event as “Benign.” **NOTE:** PREVIOUS_SIMILAR_CASES score is only calculated for copy number and Seq Var events so one must first check to make sure that the event is a copy number event before checking for the PREVIOUS_SIMILAR_CASES score or else an error will occur.

Correct usage: CASE {ANY_CN_EVENT_KIND AND SCORE(PREVIOUS_SIMILAR_CASES("Benign", 0.95)) >= 6 } { CLASSIFY("Benign") }

NOTE: If the database contains events that have previously been classified as “Benign” that are like the event in question by at least 95% and if there are at least six such similar past cases, then classify the current event as “Benign.”

Incorrect usage: CASE {SCORE(PREVIOUS_SIMILAR_CASES("Benign", 0.95)) >= 6 } { CLASSIFY("Benign") }

In the example above, a check is not made for a copy number event and therefore the statement itself is not valid and will cause an error.

SCORE(DGV)

SCORE(DGV) looks at how similar the event in question is to events in the Database of Genomic Variants (DGV). As there are likely multiple reports of polymorphisms reported in the same region, the function looks at all of these and first identifies the most similar report in terms of direction (e.g., gain or loss) as well as Similarity as defined earlier. It then takes the Similarity score and scales based on the number of samples reported in the publication. For example, if evaluating the DGV Score for a loss event, the system might identify a publication that has reported one case with 100% Similarity. In that case, the Similarity score would be approximately 65%. However, if there had been fifty such cases reported, the Similarity score would be greater than 98%.

NOTE: The DGV score is only calculated for copy number events so one must first check to make sure that the event is a copy number event before checking for the DGV score or else an error will occur.

Example usage:

Correct usage: CASE {ANY_CN_EVENT_KIND AND SCORE(DGV) > .85} {CLASSIFY("Likely Benign")}

Incorrect usage: CASE {SCORE(DGV) > .85} {CLASSIFY("Likely Benign")}

NOTE: In the example above, a check is not made for a copy number event and therefore the statement, by itself, is not valid and will return an error. The check for a copy number event must be performed each time a case statement evaluates the DGV score:

```
CASE {ANY_CN_EVENT_KIND AND SCORE(DGV) > .95} {CLASSIFY("Benign")}
CASE {ANY_CN_EVENT_KIND AND SCORE(DGV) > .85} {CLASSIFY("Likely Benign")}
```

SCORE(SAP)

The phenotypes associated with all the genes in an aberrant region are compared with the list of sample phenotypes using a statistical measure to arrive at the SAP score.

The significance of a gene based on the sample phenotype considering all phenotypes (including super classes) associated with a gene and all phenotypes associated with a sample is determined. Phenotype weights are also

considered to generate the score. A significance value is calculated using the Fisher's Exact Test and is averaged across multiple genes overlapping an event resulting in the SAP score.

SCORE(EVIDENCE)

SCORE(EVIDENCE) looks at the Evidence Score for events that have this value. The Evidence Score is calculated if a sample has HPO terms associated with it. The Evidence Score for an event is the number of genes in the region with a matching phenotype.

Example usage:

```
CASE {SCORE(EVIDENCE) >= 5} { ... }
```

NOTE: If the event contains five or more genes that match a sample phenotype, then...

```
CASE {(OVERLAP("OMIM Morbid Phenotypes Dominant") > 0) AND (SCORE(EVIDENCE) > 1)} {CLASSIFY("Likely Pathogenic")}
```

NOTE: If the event involves a dominant **OMIM Morbid** gene matching HPO terms associated with the case, then classify as "Likely Pathogenic."

Syntax for Structural Variants (SV)

EVALUATING THE TYPE OF STRUCTURAL VARIANT EVENT

The following functions assess the type of structural variant event and looks for any structural variant event:

ANY_SV_EVENT_KIND

Example usage:

```
IF { ANY_SV_EVENT_KIND } THEN { CLASSIFY("Benign") }
CASE { ANY_SV_EVENT_KIND } { CLASSIFY("Benign") }
```

SV_EVENT_KIND IS

Checks for the type of structural variant event.

Valid values for SV_EVENT_KIND:

- Insertion
- Deletion
- Inverted_duplication
- Tandem_duplication
- Inversion
- Inversion_breakpoint
- Interchr_translocation
- Intrachr_fusion

Example usage:

```
CASE { SV_EVENT_KIND IS Intrachr_fusion } {CLASSIFY("Benign") }
CASE { SV_EVENT_KIND IS Tandem_duplication } {CLASSIFY("VUS") }
CASE { SV_EVENT_KIND IS Interchr_translocation } {CLASSIFY("Likely Benign")}
```

COMPARING EVENTS TO REGION LISTS AND EVALUATING SIMILARITY SCORES

EVENT_OVERLAP

This function is like the EVENT_OVERLAP function described above.

Example usage:

```
CASE { (ANY_SV_EVENT_KIND) AND (OVERLAP("Imprinted Genes") > 0) } { CLASSIFY("Tier 1A") }
```

FUSION_MATCH

Tests to see if any of the fusion breakends (including confidence intervals) overlaps with fusion regions in a specified region list. The input to FUSION_MATCH() is a text string that is the name of an annotation list/track. A fusion is identified by two rows in the annotation list. For example, the following entries identify a MLLT3 and KMT2A translocation.

```
chr9      20341669 20622499 t(9;11)_MLLT3::KMT2A_translocation
chr11     118436492      118526832      t(9;11)_MLLT3::KMT2A_translocation
```

The FUSION_MATCH() function will then determine if the fusion breakends of an SV event in question overlap with any fusion in the region list. The function returns “true” indicating both breakends of the SV event. The function returns false if only one breakend or none of the breakends of the SV event overlap the fusions defined in the region list. This function is applicable to inter- chr translocation and intra- chr translocation and hence is applied to these functions as shown in the example below.

Example usage:

```
CASE {(SV_EVENT_KIND IS Interchr_translocation)}
{
  CASE {FUSION_MATCH("AML Translocation Interchr", 3.Kb)} {CLASSIFY("Tier 1A-Review")}
  CASE {FUSION_MATCH("Pan Heme Translocation Interchr", 3.Kb)} {CLASSIFY("Pan-Heme Overlap")}
}

CASE {(SV_EVENT_KIND IS Intrachr_fusion)}
{
  CASE {FUSION_MATCH("AML Translocation Intrachr", 3.Kb)} {CLASSIFY("Tier 1A-Review")}
  CASE {FUSION_MATCH("Pan Heme Translocation Intrachr", 3.Kb)} {CLASSIFY("Pan-Heme Overlap")}
}
```

Classifying OGM SVs

The following functions classify SV events by critical characteristics (frequency in the OGM database, SV quality and Length) to preclude these variants by being given a tiered classification and aid with filtering.

SV_PERCENTAGE_IN_OGM_CONTROLS

Filters SVs by the presence in the OGM database.

Usage:

```
CASE { ANY_SV_EVENT_KIND } {
  IF { SV_PERCENTAGE_IN_OGM_CONTROLS > 10 } THEN { CLASSIFY("Benign") }
}
```

SV_LENGTH

Filters SV by length. Best applied by specifying the SV type rather than using “ANY_SV_EVENT_KIND” .

Usage:

```
CASE { SV_EVENT_KIND IS Insertion } {
```

```

    IF { SV_LENGTH < 40.Kb } THEN { CLASSIFY("Small Variant")}
  }
  CASE { SV_EVENT_KIND IS Inverted_duplication } {
    IF { SV_LENGTH < 40. kb } THEN { CLASSIFY("Small variant")}
  }

```

SV_VAF

This function facilitates classification based on the Variant Allele Frequency (VAF) score for the selected variant.

Usage:

```

CASE { ANY_SV_EVENT_KIND } {
  IF{ SV_VAF > 0.50 } THEN { CLASSIFY("Benign") }
}

```

SV_MOLECULE_COUNT

This function facilitates classification based on the molecule count for the selected variant.

Usage:

```

CASE { ANY_SV_EVENT_KIND } {
  IF{ SV_MOLECULE_COUNT > 400 } THEN { CLASSIFY("Benign") }
}

```

SV_QUALITY

Filters SV variants by SV Quality. Quality values are imported from the VCF and are presented on a Phred scale with Bionano Solve 3.8.2. It is recommended to apply a score of Q20 for all SV types, which is equivalent to 99% accuracy in the call of the variant.

Usage:

```

CASE { ANY_SV_EVENT_KIND } {
  IF { SV_QUALITY < 20 } THEN { CLASSIFY("LowQuality") }
}

```

Evaluating Sample Attributes

ATTRIBUTE()

This function evaluates values of specified factor names (attributes). Each sample can have several attributes, e.g., Gender, Tumor Type, Age, etc.

Allowed operators:

When comparing using numbers,

- > Greater than
- < Less than
- >= Greater than or equal to
- != Not equal to
- <= Less than or equal to
- == Equal to (**NOTE:** must use two equal signs in succession)

When the value being evaluated is a text string, the following operators are used:

- == Equals (returns a “true” or “false”)
- != Does not equal

Example usage:

```
CASE {ATTRIBUTE("Gender") == "male" } {...}
CASE {ATTRIBUTE("Age") < 30} {...}
```

Example Decision Tree Script

Below is a sample decision tree script using CASE statements for pre-classification of constitutional analysis. This script can be used as a template and further customized with in-house tracks from user past cases.

```
/*
Using ACMG Guidelines for CNV calling of events not previously classified as Benign, Likely Benign, or Artifact more
than 4 times
Classifies CNV events as:
Benign
Likely Benign
VUS
Likely Pathogenic
Pathogenic
Artifact
Exclude

Classifies AOH events as
*/
val Sim_threshold = 0.85

CASE {ANY_CN_EVENT_KIND} {

CASE {(SCORE(PREVIOUS_SIMILAR_CASES("Pathogenic", 0.85)) == 0) AND
```

```
(SCORE(PREVIOUS_SIMILAR_CASES("Likely Pathogenic", 0.85)) == 0) AND
(SCORE(PREVIOUS_SIMILAR_CASES("VUS", 0.85)) == 0)} {
CASE {SCORE(PREVIOUS_SIMILAR_CASES("Artifact", 0.95)) >=4} {CLASSIFY("Artifact")}
CASE {SCORE(PREVIOUS_SIMILAR_CASES("Benign", 0.95)) >=4} {CLASSIFY("Benign")}
CASE {SCORE(PREVIOUS_SIMILAR_CASES("Likely Benign", 0.95)) >=4} {CLASSIFY("Likely Benign")}
}
CASE {SCORE(PREVIOUS_SIMILAR_CASES("Pathogenic", 0.95)) >=1} {CLASSIFY("Pathogenic")}
CASE {SCORE(PREVIOUS_SIMILAR_CASES("Likely Pathogenic", 0.95)) >=4} {CLASSIFY("Likely Pathogenic")}
}

CASE {OVERLAP("RefSeq Genes") == 0} {CLASSIFY("Exclude")}

CASE {ANY_CN_EVENT_KIND AND SCORE(DGV)>0.88} {CLASSIFY("Likely Benign")} // Classify likely benign based on DGV
cases/score and absence of COE Paper genes as Likely benign No Reporting

CASE {ANY_CN_EVENT_KIND} {

CASE{SCORE(ACMG_CNV) >= 0.99} {CLASSIFY("Pathogenic")}
CASE{SCORE(ACMG_CNV) >= 0.9} {CLASSIFY("Likely Pathogenic")}
CASE{SCORE(ACMG_CNV) >= -0.89} {CLASSIFY("VUS")}
CASE{SCORE(ACMG_CNV) >= -0.98} {CLASSIFY("Likely Benign")}
CASE{SCORE(ACMG_CNV) < -0.98} {CLASSIFY("Benign")}
}
}
```

Tiered Variant Decision Tree for Hematological Malignancies

The latest release of VIA includes reference files with the annotation and tracks update for guideline targets, which are region coordinates that are derived from various global societies specific to Acute Myeloid Leukemia (AML), Myelodysplastic Syndrome (MDS), and an aggregation of all hematological malignancy guidelines (pan-heme). These resources, along with custom region files, can be used for a customizable decision tree for CNV, AOH, and SV event pre-classification. New installations of VIA include a standard decision tree for these three disease states that is constructed by Bionano to facilitate the identification of relevant events for researchers.

Schematic for Disease-specific Decision Tree

A preconfigured decision tree is provided for customization by the site administrators, as necessary. Decision trees are constructed to preclassify SV and CNV events based on overlap to Heme-guideline target lists curated by Bionano. Variant pre-classification in VIA with the rules-based decision tree is serial, meaning the order of syntax is important as an event proceeds through the decision tree linearly until it is classified, or ultimately left

unclassified. The classification schema outlined in **Figure 27** is an example and fully configurable to leverage differing classification states, resources, or tiering strategy to achieve site-specific objectives.

		SV Event Types							CN Event Types				LOH					
		Del	Ins	Inv	Inter	Intra	Tan Dup	Inv Dup	Loss	Null	Gain	Amp	AOH					
Disease Specific	Guideline Target Files	Tier 1A	Tier 1A	Tier 1A	Tier 1A - Review	Tier 1A - Review	Tier 1A	Tier 1A	Artifact Tier 1A	Artifact Tier 1A	Artifact	Artifact						
	CNV Solve Mask																	
	Deletion- Large																	
	Deletion- Small																	
	Duplication- Large																	
	Duplication- Small																	
	Inversion																	
	Interchrom Trans																	
	Intrachrom Trans																	
	Rearrangements - Large																	
	Rearrangements -Small	Tier 1A	Tier 1A	Tier 1A	Tier 1A	Tier 1A	Tier 1A	Tier 1A	Tier 1A	Tier 1A	Tier 1A							
	Monosomy																	
	Trisomy																	
Pan Heme	Deletion- Large	Pan-Heme Overlap	Pan-Heme Overlap	Pan-Heme Overlap	Pan-Heme Overlap	Pan-Heme Overlap	Pan-Heme Overlap	Pan-Heme Overlap	Pan-Heme Overlap	Pan-Heme Overlap								
	Deletion- Small																	
	Duplication- Large																	
	Duplication- Small																	
	Inversion																	
	Interchrom Trans																	
	Intrachrom Trans																	
	Rearrangements - Large																	
	Rearrangements -Small	Pan-Heme Overlap							Pan-Heme Overlap	Pan-Heme Overlap	Pan-Heme Overlap	Pan-Heme Overlap	Pan-Heme Overlap	Pan-Heme Overlap	Pan-Heme Overlap	Pan-Heme Overlap	Pan-Heme Overlap	
	Monosomy																	
	Trisomy																	
	CIVIC	CIVIC Genes							CIVIC Gene	CIVIC Gene	CIVIC Gene	CIVIC Gene	CIVIC Gene	CIVIC Gene	CIVIC Gene	CIVIC Gene	CIVIC Gene	CIVIC Gene
Copy Number	Prev to Sim Tier 1B								Tier 1B	Tier 1B	Tier 1B	Tier 1B						
	Prev to Sim Tier 2								Tier 2	Tier 2	Tier 2	Tier 2						
	Prev to Sim Tier 3								Tier 3	Tier 3	Tier 3	Tier 3						
	Prev to Sim Tier 4								Tier 4	Tier 4	Tier 4	Tier 4						
	Prev to Sim Artifact								Artifact	Artifact	Artifact	Artifact						
	DGV Similar >88%								Tier 4	Tier 4	Tier 4	Tier 4						
	Unclassified																	

Figure 27. Legend: Del – Deletion, Ins = Insertion, Inv = Inversion, Inter = Inter-chromosomal translocation, Intra = Intrachromosomal fusion, Tan Dup = Tandem Duplication, Inv Dup = Inverted Duplication, Null = Homozygous Loss, Amp = Amplification, LOH = Copy Neutral Loss of Heterozygosity

Decision Tree Script For Disease Specific Tiering

Example syntax for decision tree to pre-classify CNV and SV variants obtained with OGM data with Bionano guideline files specific for AML as illustrated in the Decision Tree schematic above.

```
//
Event Classifications:
Tier 1A
Tier 1B
Tier 2
Tier 3
Tier 4
```

Tier 1A-Review
Pan-Heme Overlap
CIVIC Gene
Artifact

```

*/
// CNV Masking overlap Artifacts
CASE {(ANY_CN_EVENT_KIND)}
{
  CASE {EVENT_OVERLAP("Solve CNV Mask") > 0.45} {CLASSIFY("Artifact")}
}

// Gene Fusions disease-specific
CASE {(SV_EVENT_KIND IS Interchr_translocation)}
{
  CASE {FUSION_MATCH("AML Translocation Interchr", 3.Kb)} {CLASSIFY("Tier 1A-Review")}
}

CASE {(SV_EVENT_KIND IS Intrachr_fusion)}
{
  CASE {FUSION_MATCH("AML Translocation Intrachr", 3.Kb)} {CLASSIFY("Tier 1A-Review")}
}

// Deletions disease-specific
CASE {(SV_EVENT_KIND IS Deletion)}
{
  CASE {EVENT_OVERLAP("AML Deletion Small", 3.Kb) > 0} {CLASSIFY("Tier 1A")}
}

CASE {(CN_EVENT_KIND IS CN_LOSS) OR (CN_EVENT_KIND IS HOMOZYGOUS_COPY_LOSS)}
{
  CASE {SIMILARITY("AML Deletion Small") > 0.9} {CLASSIFY("Tier 1A")}
  CASE {SIMILARITY("AML Deletion Large") > 0.9} {CLASSIFY("Tier 1A")}
}

// **Duplication disease-specific
CASE {(SV_EVENT_KIND IS Tandem_duplication) OR (SV_EVENT_KIND IS Inverted_duplication) OR (SV_EVENT_KIND IS
  Insertion)}
{
  CASE {EVENT_OVERLAP("AML Duplication Small", 3.Kb) > 0} {CLASSIFY("Tier 1A")}
}

```

```

    }

CASE {(CN_EVENT_KIND IS CN_GAIN) OR (CN_EVENT_KIND IS HIGH_COPY_GAIN)}
{
    CASE {SIMILARITY("AML Duplication Small") > 0.9} {CLASSIFY("Tier 1A")}
}

// **Monosomy disease-specific
CASE {(CN_EVENT_KIND IS CN_LOSS)}
{
    CASE {SIMILARITY("AML Monosomy") > 0.9} {CLASSIFY("Tier 1A")}
}

// **Trisomy disease-specific
CASE {(CN_EVENT_KIND IS CN_GAIN)}
{
    CASE {SIMILARITY("AML Trisomy") > 0.9} {CLASSIFY("Tier 1A")}
}

// **Rearrangements disease-specific
CASE {(ANY_SV_EVENT_KIND)}
{
    CASE {EVENT_OVERLAP("AML Rearrangements Small", 3.Kb) > 0} {CLASSIFY("Tier 1A")}
}

CASE {(ANY_CN_EVENT_KIND)}
{
    CASE {SIMILARITY("AML Rearrangements Small") > 0.9} {CLASSIFY("Tier 1A")}
}

// Gene Fusions Pan-Heme
CASE {(SV_EVENT_KIND IS Interchr_translocation)}
{
    CASE {FUSION_MATCH("Pan Heme Translocation Interchr", 3.Kb)} {CLASSIFY("Pan-Heme Overlap")}
}

```

```

CASE {(SV_EVENT_KIND IS Intrachr_fusion)}

{

CASE {FUSION_MATCH("Pan Heme Translocation Intrachr", 3.Kb)} {CLASSIFY("Pan-Heme Overlap")}

}


// Deletions Pan-Heme

CASE {(SV_EVENT_KIND IS Deletion)}

{

CASE {EVENT_OVERLAP("Pan Heme Deletion Small", 3.Kb) > 0} {CLASSIFY("Pan-Heme Overlap")}

}


CASE {(CN_EVENT_KIND IS CN_LOSS) OR (CN_EVENT_KIND IS HOMOZYGOUS_COPY_LOSS)}

{

CASE {SIMILARITY("Pan Heme Deletion Small") > 0.9} {CLASSIFY("Pan-Heme Overlap")}

CASE {SIMILARITY("Pan Heme Deletion Large") > 0.9} {CLASSIFY("Pan-Heme Overlap")}

}


// **Duplication Pan-Heme

CASE {(SV_EVENT_KIND IS Tandem_duplication) OR (SV_EVENT_KIND IS Inverted_duplication) OR (SV_EVENT_KIND IS
Insertion)}

{

CASE {EVENT_OVERLAP("Pan Heme Duplication Small", 3.Kb) > 0} {CLASSIFY("Pan-Heme Overlap")}

}


CASE {(CN_EVENT_KIND IS CN_GAIN) OR (CN_EVENT_KIND IS HIGH_COPY_GAIN)}

{

CASE {EVENT_OVERLAP("Pan Heme Duplication Small") > 0} {CLASSIFY("Pan-Heme Overlap")}

CASE {EVENT_OVERLAP("Pan Heme Duplication Large") > 0} {CLASSIFY("Pan-Heme Overlap")}

}


// **Monosomy Pan-Heme

CASE {(CN_EVENT_KIND IS CN_LOSS)}

{

CASE {SIMILARITY("Pan Heme Monosomy") > 0.9} {CLASSIFY("Pan-Heme Overlap")}

}

```

```
// **Trisomy Pan-Heme
CASE {(CN_EVENT_KIND IS CN_GAIN)}
{
    CASE {SIMILARITY("Pan Heme Trisomy") > 0.9} {CLASSIFY("Pan-Heme Overlap")}
}

// **Rearrangements Pan-Heme
CASE {(ANY_SV_EVENT_KIND)}
{
    CASE {EVENT_OVERLAP("Pan Heme Rearrangements Small", 3.Kb) > 0} {CLASSIFY("Pan-Heme Overlap")}
}

CASE {(ANY_CN_EVENT_KIND)}
{
    CASE {SIMILARITY("Pan Heme Rearrangements Small") > 0.9} {CLASSIFY("Pan-Heme Overlap")}
    CASE {SIMILARITY("Pan Heme Rearrangements Large") > 0.9} {CLASSIFY("Pan-Heme Overlap")}
}

//Cancer Genes CNV
CASE {(ANY_CN_EVENT_KIND)}
{
    CASE {EVENT_OVERLAP("CIViC Genes")>0} {CLASSIFY("CIViC Gene")}
}

//Cancer Genes SV
CASE {(ANY_SV_EVENT_KIND)}
{
    CASE {EVENT_OVERLAP("CIViC Genes", 3.Kb)>0} {CLASSIFY("CIViC Gene")}
}

//Case History
CASE {(ANY_CN_EVENT_KIND)}
{
    CASE {SCORE(PREVIOUS_SIMILAR_CASES("Tier 1B", 0.95)) >=1} {CLASSIFY("Tier 1B")}
}
```

```

CASE {SCORE(PREVIOUS_SIMILAR_CASES("Tier 2", 0.95)) >=1} {CLASSIFY("Tier 2")}
}

CASE {(ANY_CN_EVENT_KIND) AND
(SCORE(PREVIOUS_SIMILAR_CASES("Tier 1A", 0.85)) == 0) AND
(SCORE(PREVIOUS_SIMILAR_CASES("Tier 1B", 0.85)) == 0) AND
(SCORE(PREVIOUS_SIMILAR_CASES("Tier 2", 0.85)) == 0)}
{
CASE {SCORE(PREVIOUS_SIMILAR_CASES("Tier 3", 0.95)) >=4} {CLASSIFY("Tier 3")}
CASE {SCORE(PREVIOUS_SIMILAR_CASES("Tier 4", 0.95)) >=4} {CLASSIFY("Tier 4")}
CASE {SCORE(PREVIOUS_SIMILAR_CASES("Artifact", 0.95)) >=4} {CLASSIFY("Artifact")}
}

// DGV
CASE {ANY_CN_EVENT_KIND AND SCORE(DGV) > 0.88} {CLASSIFY("Tier 4")}

```

OPTIONAL DECISION TREE SCRIPT ADDITION FOR VIA 7.1 +

The following text string can be added to the decision tree to utilize the new SV autoclassification functionality. The local administrator must create new **Event Classification** states for the corresponding designations.

```

// Labeling low qualitySVs
}

//SV Quality
CASE {ANY_SV_EVENT_KIND AND SV_QUALITY < 20} { CLASSIFY("Low Quality SV")}

//SV Control database for heme - 0%, use 1 for constitutional
CASE { ANY_SV_EVENT_KIND AND SV_PERCENTAGE_IN_OGM_CONTROLS > 0 } { CLASSIFY("In Control
Database")}

//SV length - small SVs
CASE {(SV_EVENT_KIND IS Deletion AND SV_LENGTH < 2.Kb) OR (SV_EVENT_KIND IS Tandem_duplication AND
SV_LENGTH < 2.Kb) OR (SV_EVENT_KIND IS Inverted_duplication AND SV_LENGTH < 2.Kb) OR
(SV_EVENT_KIND IS Insertion AND SV_LENGTH < 2.Kb) OR (SV_EVENT_KIND IS Inversion AND SV_LENGTH <
2.Kb)} { CLASSIFY("SV < 2 kbp")}

```

Appendix

AOH performance settings

```
#Step = BAF from OGM.BAM
# Type = BAF from OGM BAM Processing
# Reject labels with coverage less than = 20
# Reject reads with MAPQ less than = 30
# OGM cluster file for BAF = GmmMaf05_220623.tsv.gz
#Step = CN from BAM MultiScale
# Type = Multiscale Reference BAM Processing
#Step = Read Depth Track
# Type = BAM Read Depth Track Analysis
#Processing Type Parameters:
#Step = Recenter Probes
# Allow manual centering = true
# Type = Median
#Step = Analysis
# Type = SNP-FASST3 Segmentation
# Max Contiguous Probe Spacing (Kbp) = 1000.0
# Min number of probes per segment = 3
# CN Significance Threshold = 5.0E-6
# SNP Significance Threshold = 0.005
# Segment Boundaries = Midpoint between adjacent probes
# Amplification (4+:2) = 0.6
# Gain (3:2) = 0.2
# Loss (1:2) = -0.2
# Homozygous Loss (0:2) = -1.0
# M vs M, X/Y Loss (0:1) = -1.0
# M vs F, X Loss (0:2) = -1.0
# F vs M, X Loss (1:1) = 0.2
# M vs M, X/Y Gain (2:1) = 0.2
# M vs F, X Gain (2:2) = -0.2
# M vs M, X/Y Amplification, F vs M, X Gain (3:1) = 1.2
# F vs M, X Amplification (4+:1) = 1.7
# Homozygous Frequency Threshold = 0.85
# Homozygous Value Threshold = 0.9
# Heterozygous Imbalance Threshold = 0.4
# Minimum LOH Length (KB) = 500
# Minimum SNP Probe Density (Probes/MB) = 0.0
# Min. mosaic threshold (%) = 15
# Mosaic CN significance offset = 1.0E-12
# Mosaic SNP significance offset = 0.0
#Step = Robust Variance Sample QC Calculation
```

```
# Percent outliers to remove = 3.0
#Step = Mosaic Analysis
# Type = Aberrant Cell Fraction
# Mosaic Labeling = true
# Label anything lower than (%) = 80.0
#Array Decision Tree = None
#String Filter:
#CN Filters:
#Allelic Filters:
# With Copy Number Events:
# Remove AOH regions that are also copy number loss = yes
# Remove allelic imbalance regions that are also copy number events = yes
```

Technical Assistance

For technical assistance, contact Bionano Technical Support.

You can retrieve documentation on Bionano products, SDS's, certificates of analysis, frequently asked questions, and other related documents from the Support website or by request through e-mail and telephone.

TYPE	CONTACT
Email	support@bionano.com
Phone	Hours of Operation: Monday through Friday, 9:00 a.m. to 5:00 p.m., PST US: +1 (858) 888-7663 Monday through Friday, 9:00 a.m. to 5:00 p.m., CET UK: +44 115 654 8660 France: +33 5 37 10 00 77 Belgium: +32 10 39 71 00
Website	www.bionano.com/support
Address	Bionano, Inc. 9540 Towne Centre Drive, Suite 100 San Diego, CA 92121

Legal Notice

For Research Use Only. Not for use in diagnostic procedures.

This material is protected by United States Copyright Law and International Treaties. Unauthorized use of this material is prohibited. No part of the publication may be copied, reproduced, distributed, translated, reverse-engineered or transmitted in any form or by any media, or by any means, whether now known or unknown, without the express prior permission in writing from Bionano Genomics. Copying, under the law, includes translating into another language or format. The technical data contained herein is intended for ultimate destinations permitted by U.S. law. Diversion contrary to U. S. law prohibited. This publication represents the latest information available at the time of release. Due to continuous efforts to improve the product, technical changes may occur that are not reflected in this document. Bionano Genomics reserves the right to make changes to specifications and other information contained in this publication at any time and without prior notice. Please contact Bionano Genomics Customer Support for the latest information.

BIONANO GENOMICS DISCLAIMS ALL WARRANTIES WITH RESPECT TO THIS DOCUMENT, EXPRESSED OR IMPLIED, INCLUDING BUT NOT LIMITED TO THOSE OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE. TO THE FULLEST EXTENT ALLOWED BY LAW, IN NO EVENT SHALL BIONANO GENOMICS BE LIABLE, WHETHER IN CONTRACT, TORT, WARRANTY, OR UNDER ANY STATUTE OR ON ANY OTHER BASIS FOR SPECIAL, INCIDENTAL, INDIRECT, PUNITIVE, MULTIPLE OR CONSEQUENTIAL DAMAGES IN CONNECTION WITH OR ARISING FROM THIS DOCUMENT, INCLUDING BUT NOT LIMITED TO THE USE THEREOF, WHETHER OR NOT FORESEEABLE AND WHETHER OR NOT BIONANO GENOMICS IS ADVISED OF THE POSSIBILITY OF SUCH DAMAGES.

Patents

Products of Bionano Genomics® may be covered by one or more U.S. or foreign patents.

Trademarks

The Bionano logo and names of Bionano products or services are registered trademarks or trademarks owned by Bionano Genomics, Inc. ("Bionano") in the United States and certain other countries.

Bionano™, Bionano Genomics®, Bionano Access™, and VIA™ are trademarks of Bionano Genomics, Inc. All other trademarks are the sole property of their respective owners.

No license to use any trademarks of Bionano is given or implied. Users are not permitted to use these trademarks without the prior written consent of Bionano. The use of these trademarks or any other materials, except as permitted herein, is expressly prohibited and may be in violation of federal or other applicable laws.

© Copyright 2025 Bionano Genomics, Inc. All rights reserved.